



## Video microscope

OIV 355

TOIV 355-A



PROFESSIONAL MEASURING

### Operating instructions Video microscope OIV 355

Version 1.0  
2026-01  
GB  
OIV-35-BA-d-2610

**KERN Optics OIV 355**

**Video microscope**

**Operating instructions Video microscope OIV 355**

Version 1.0 2026-01

Table of contents:

<b>1</b>	<b>Technical data</b> .....	<b>3</b>
<b>2</b>	<b>Declaration of conformity</b> .....	<b>4</b>
<b>3</b>	<b>Overview of the device</b> .....	<b>5</b>
3.1	Scope of delivery .....	5
3.2	Nomenclature .....	6
<b>4</b>	<b>Basic information (general)</b> .....	<b>8</b>
4.1	General information on warnings .....	8
4.2	Intended use .....	9
4.3	Improper use .....	9
4.4	Warranty .....	9
<b>5</b>	<b>Basic warnings and safety instructions</b> .....	<b>10</b>
5.1	Observe the notes in the operating instructions .....	10
5.2	Staff training .....	10
5.3	Security .....	10
<b>6</b>	<b>Transport and storage</b> .....	<b>12</b>
6.1	Note .....	12
6.2	Transport .....	12
6.3	Storage .....	12
6.4	Packaging/return transport .....	12
<b>7</b>	<b>Unboxing and commissioning</b> .....	<b>13</b>
7.1	Unboxing .....	13
<b>8</b>	<b>Assembly</b> .....	<b>14</b>
8.1	Assembling the transmitted light stand.....	14
8.2	Mounting the microscope head on the stand .....	14
8.3	Connecting to the power supply .....	15
<b>9</b>	<b>Operation and function of the microscope</b> .....	<b>16</b>
9.1	Switching on the device.....	16
9.2	Switching on the lights.....	16
9.3	Setting the illuminance (incident & transmitted light) .....	16
9.4	Adjusting magnification and focus.....	17
9.5	Connections and interfaces .....	17
9.6	Inserting the MicroSD memory card.....	18
<b>10</b>	<b>Operation and function of the microscope</b> .....	<b>19</b>
10.1	Camera menu.....	19
10.2	Exposure and white balance .....	19
10.3	Adjust parameters.....	20
10.4	Control elements .....	21
10.5	Measurements .....	23
10.6	Measuring tools .....	23
10.7	Calibration setting.....	27
10.8	Measurement attributes.....	29

10.9	Scale ruler.....	30
10.10	Further settings.....	30
10.11	Languages.....	30
10.12	Storage space .....	31
10.13	Resetting the data to the factory settings .....	31
10.14	Set date / time.....	32
10.15	Switch between two displays.....	32
10.16	Move and hide toolbar .....	32
<b>11</b>	<b>Toolbar .....</b>	<b>33</b>
11.1	Camera.....	33
11.2	Camera timer.....	33
11.3	Video recording .....	35
11.4	Gallery .....	35
11.5	Rename image and video files .....	36
11.6	Delete image and video files .....	37
11.7	Open image files.....	37
11.8	Open video files.....	37
11.9	File details .....	37
<b>12</b>	<b>Optional accessories .....</b>	<b>38</b>
<b>13</b>	<b>Troubleshooting .....</b>	<b>39</b>
<b>14</b>	<b>Service.....</b>	<b>40</b>
<b>15</b>	<b>Power supply .....</b>	<b>41</b>
15.1	Opening video files.....	41
<b>16</b>	<b>Maintenance, servicing and disposal.....</b>	<b>42</b>
16.1	Maintenance and cleaning .....	42
16.2	Repair .....	42
16.3	Disposal of the appliance and accessories .....	42
<b>17</b>	<b>Battery Act .....</b>	<b>44</b>
<b>18</b>	<b>Further information .....</b>	<b>45</b>

## 1 Technical data

odel KERN	OIV 355
Item number/type	TOIV 355-A
Optical system	Greenough System
Lighting	LED ring light 3 W, LED transmitted light 5 W
Dimmable lighting	Yes
Stand / height / base plate	Mechanical / 330 mm / 320 x 260 mm
Zoom Optical / Digital	0.7x – 4.5x / 17x – 110x
Magnification ratio	6.4 : 1
Field of vision Ø	17.6 - 1.47 mm
Screen tilt angle	48°
Screen size / resolution	11.6 inches / 1920 x 1080 Full HD
Camera sensor / resolution	SONY CMOS 1/1.8 Zoll / 5 MP
Refresh rate	60 fps at 1920x x 1080
Recording	60 fps at 1920x x 1080
Interfaces	HDMI, USB 2.0 (2x)
Storage of data	Snapshot, recording, USB
Camera Parameters	Exposure, white balance, gain, gamma, brightness, contrast, hue, saturation, greyscale
Measurements	Point, line, parallel lines, dashed lines, angle, diameter circle, radius circle, three-point circle, distance from circle to line, concentric circles, double circle, vertical line, any line, polygon, arc, marking, screenshot, table data export, calibration, parameter settings
Further functions	Horizontal and vertical mirroring, freezing
Languages	German, English, French, Italian, Russian, Japanese, Chinese
Net / gross weight	4.5 kg / 6 kg
Dimensions packaging	420 x 390 x 250 mm

## **2 Declaration of conformity**

The current EC/EU Declaration of Conformity can be found online

### **3 Overview of the device**

#### **3.1 Scope of delivery**

- Video microscope
- Power supply unit 100-240 V, 50/60 Hz, output 12 V DC, 2000 mA
- USB wireless mouse
- USB memory stick (32 GB)
- Calibration slide (multi)
- Stand insert frosted glass, 95 mm
- Stand insert plastic, black/white, 95 mm
- Operating instructions

3.2 Nomenclature





10



11

12

Position	Description
1	Microscope stand
2	LED ring light
3	Microscope head with camera
4	Screen
5	Zoom adjustment wheel with click-stop function
6	Focussing wheel
7	Dimmer LED ring light
8	Dimmer LED transmitted light unit
9	Transmitted light unit with stand insert & object holder
10	Power supply input and output
11	Interfaces
12	Memory card slot





## 4 Basic information (general)


### 4.1 General information on warnings

Warnings are used in these operating instructions to warn you of possible personal injury or damage to property in certain situations.

Signal words	Description
<b>DANGER</b>	Failure to observe the instructions will lead directly to serious injury, permanent impairment (e.g. loss of a limb) or death of the user or third parties
<b>WARNING</b>	Failure to observe the instructions may result in serious injury, permanent impairment (e.g. loss of a limb) or death of the user or third parties
<b>CAUTION</b>	Failure to observe the instructions may result in minor injuries or temporary damage to the user or third parties (e.g. minor cuts)
<b>NOTE</b>	Failure to observe the instructions may result in damage to property

#### Symbols and warnings:

Symbol	Meaning
<b>Warning signs</b>	Warning signs warn you of hazards that may result in personal injury. The symbol indicates the type of hazard.
	Indicates general hazards or a danger point
	Warning of electrical voltage

Symbol	Meaning
<b>Command sign</b>	Mandatory signs prescribe measures that you must take to avoid personal injury or damage to property. The symbol indicates the necessary actions or objects to prevent damage.
	Indicates a prescribed action

## **4.2 Intended use**

The intended use of the OIV 355 video microscope is the non-contact, optical magnification and digital visualisation of objects or surfaces for visual inspection, documentation or analysis. It is typically used in laboratories, quality assurance or research to visualise details that cannot be seen with the naked eye. It is operated via a camera unit that transmits the image to a monitor.

## **4.3 Improper use**

The OIV 355 video microscope is not to be used for medical purposes. Do not use the device in potentially explosive atmospheres or for measurements in liquids or on live parts. Unauthorised structural changes, additions or conversions to the device are prohibited.

## **4.4 Warranty**

The warranty is void if

- Non-compliance with our specifications in the operating instructions
- Use outside the described applications
- Modification or opening of the device
- Mechanical damage and damage caused by media, liquids, natural wear and tear
- Improper installation, incorrect assembly or electrical installation

## 5 Basic warnings and safety instructions

### 5.1 Observe the notes in the operating instructions




Read the operating instructions carefully before commissioning/using the device, even if you already have experience with KERN devices. Always keep the manual in the immediate vicinity of the device. If you lose these instructions, you can find them in our download area at [www.kern-sohn.com](http://www.kern-sohn.com)

### 5.2 Staff training

The appliance may only be used by persons who have read and understood the operating instructions, particularly the chapter on safety.

### 5.3 Security

<b>⚠️ WARNUNG</b>	
	<p>Read all safety information and instructions. Failure to comply with the safety information and instructions may result in electric shock, fire and/or serious injury. Keep all safety information and instructions for future reference.</p> <ul style="list-style-type: none"><li>• The microscope must not be modified in any way. This can lead to incorrect measurement results, safety-related defects and the destruction of the measuring device</li><li>• Do not operate the device in potentially explosive rooms or areas and do not set it up there.</li><li>• Do not operate the device in an aggressive atmosphere.</li><li>• Do not immerse the device in water. Do not allow any liquids to penetrate the inside of the device.</li><li>• The device may only be used in a dry environment and under no circumstances in rain or at a relative humidity above the operating conditions.</li><li>• Protect the device from permanent direct sunlight.</li><li>• Do not expose the device to strong vibrations.</li><li>• Do not remove any safety signs, stickers or labels from the device. Keep all safety signs, stickers and labels in a legible condition.</li><li>• Do not open the appliance.</li></ul>

### **⚠ WARNING**



**Risk of injury due to electric shock!**

- Risk of short circuit due to liquids penetrating the housing!
- Do not immerse the device and accessories in water. Ensure that no water or other liquids get into the housing.
- Work on electrical components may only be carried out by an authorised specialist company!
- Take care not to twist or kink the mains or connection cable.
- Only use the original adapter supplied.

### **⚠ WARNING**



**Risk of suffocation!**

Do not leave the packaging material lying around carelessly. It could become a dangerous toy for children.

- The appliance is not a toy and does not belong in the hands of children.
- This appliance can be dangerous if it is used improperly or not as intended by untrained persons! Therefore, observe the personnel qualifications!

### **⚠ WARNING**



**Electrostatic sensitive device!**

- The device can be destroyed by electrostatic discharge. Connectors for HF signals are particularly at risk.
- Please observe the handling instructions for electrostatically sensitive components.

### **CAUTION**

- Keep a sufficient distance from heat sources.
- Do not use the device in environments with high humidity or water mist.

### **! NOTE**

- To avoid damaging the appliance, do not expose it to extreme temperatures, extreme humidity or moisture.
- Do not use harsh cleaning agents, abrasive cleaners or solvents to clean the appliance.

## **6 Transport and storage**

### **6.1 Note**

If you store or transport the device improperly, the device may be damaged. Observe the information on transporting and storing the device.

Avoid exposing the device to direct sunlight, high or low temperatures, shocks, dust and high humidity.

The suitable temperature range is 0 - 40 °C and a relative humidity of 85 % should not be exceeded.

The device should always be placed on a firm, smooth and horizontal surface.

### **6.2 Transport**

Always use the original packaging when transporting the appliance to protect it from external influences.

### **6.3 Storage**

Observe the following storage conditions when the device is not in use:

- dry and protected from frost and heat
- Cover the microscope with the dust cover supplied
- the storage temperature corresponds to the technical data

### **6.4 Packaging/return transport**

Returns are only possible within the limits of the general terms and conditions. Keep all parts of the original packaging for any necessary return transport.

- Only use the original packaging for return transport.
- Disconnect all connected cables and loose/movable parts before despatch.
- Reattach any transport locks provided.
- Secure all parts against slipping and damage.

## **7 Unboxing and commissioning**

### **7.1 Unboxing**



In the event of a return, please follow the instructions in the "Packaging/return transport" chapter.

---

On receipt of the appliance, you should first check that no damage has occurred during transport, that the outer packaging, the housing, other parts or even the appliance itself have not been damaged. If any damage is evident, please notify KERN & SOHN GmbH immediately.

## 8 Assembly

### 8.1 Assembling the transmitted light stand




Attach the stand column to the base plate by inserting the three screws supplied (incl. spring washer, see illustration) through the plate and turning them into the underside of the column using the Allen key supplied.

Tighten the screws hand-tight.

Set up the microscope stand after successful assembly.

### 8.2 Mounting the microscope head on the stand

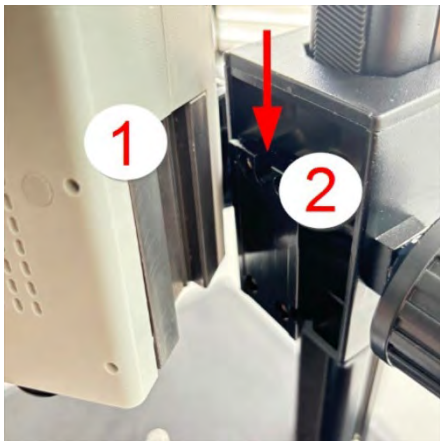

Then insert the microscope head (1) into the guide on the holder (2) of the microscope stand in the direction of the arrow.



Then tighten the two inside screws on using the socket spanner supplied.

Make sure that the socket spanner has a good grip.


Tighten the screws hand-tight.

### 8.3 Connecting to the power supply



Insert the angled power supply plug into the power input socket/input (1).

Next, take the spiral cable supplied and plug one end of it into the output socket (2).

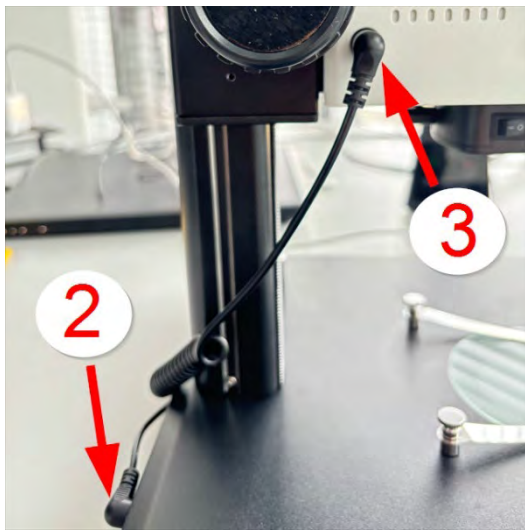


Then connect the other end of the spiral cable to the microscope head (3).

Note: Ensure that all plug connections are secure. Carry out this check again, especially after a change of location.



## 9 Operation and function of the microscope

### 9.1 Switching on the device

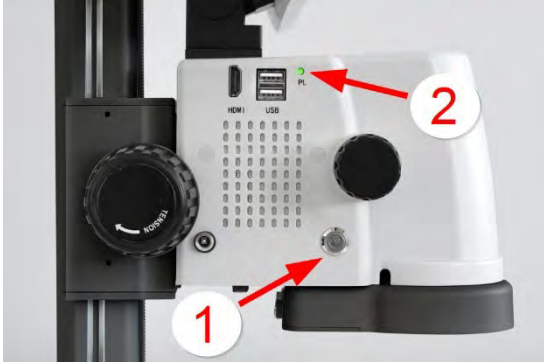

To switch on the video microscope, press the On/Off button (1) until it is backlit in blue.

The LED operating indicator (2) lights up green as soon as the microscope is switched on and switches to operating mode.

### 9.2 Switching on the lights



The on/off switch for the LED ring lighting is located on the side below the main switch.

### 9.3 Setting the illuminance (incident & transmitted light)



The control for adjusting the incident light illumination is located on the opposite side of the on/off switch. Turn the control in the Min or Max direction to set the light intensity perfectly.



The intensity of the transmitted light can be found on the top of the base plate, in the back right corner. You can use this control to adjust the brightness separately from the ring light.


## 9.4 Adjusting magnification and focus



As the OIV series are zoom microscopes, the zoom factor (0.7x-4.5x) can be changed using the zoom adjustment wheel (1). Turn the adjustment wheel to the maximum magnification. Focus the image using the focussing wheel (2). Repeat the process if necessary.

Note: The max/min parfocality is already optimally preset on delivery.

## 9.5 Connections and interfaces

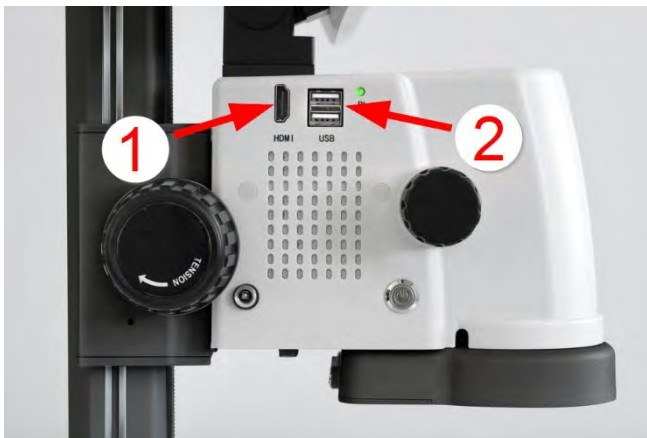

The device has various interfaces:

- (1) HDMI
- (2) USB 2.0 (2x)

These are used for different applications of the internal and external software.

### Internal software via the fixed screen:

To operate the software, a mouse must be connected via the USB port. To do this, plug the USB receiver, which is located in the battery compartment of the wireless mouse supplied, into one of the two USB slots (2). A USB stick (max. 128 GB) must be connected for data storage. Note: There is no internal storage capacity.

### External software via a PC / laptop:

The video microscope can be connected to a PC/laptop using a USB cable (type A to type A). To do this, connect the USB cable to one of the two USB slots (2) and to a USB port on the PC/laptop.

The external software (S-Viewer) is available free of charge under Downloads on our website >>> <https://www.kern-sohn.com/shop/de/downloads>

The image can be transmitted in parallel to another (larger) display device via the additional HDMI output (1).

## 9.6 Inserting the MicroSD memory card




There is a slot for a microSD memory card on the back of the device. The connection thus offers an additional option for storing data (max. 128 GB).

To do this, insert a MicroSD memory card (not included) into the card slot until you hear and feel it click into place.


## 10 Operation and function of the microscope



### 10.1 Camera menu

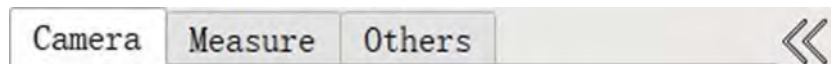

Click on the **camera** to switch to the camera page, set the relevant camera parameters and use the corresponding tools.




### 10.2 Exposure and white balance

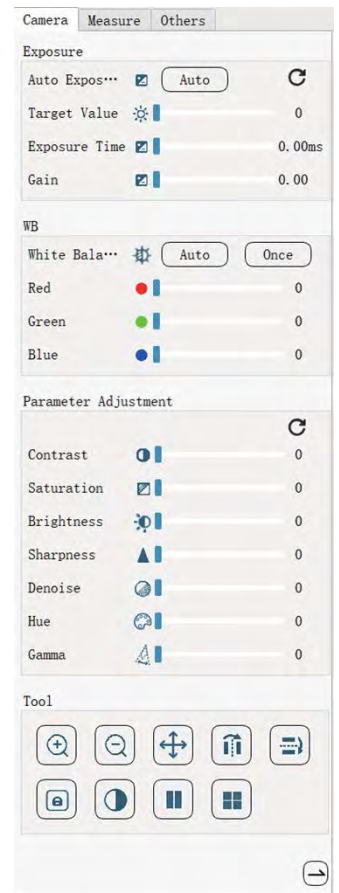

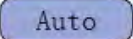
 Automatic exposure is activated. The target value can be set by dragging the  slider to the left or right.

This deactivates manual adjustment of the exposure time or the amplifier (gain)

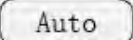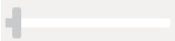


 Automatic exposure is deactivated. The target value cannot be changed  Exposure time and amplifier (gain) can be  set manually using the sliders.

 Click this button to reset the exposure parameters and activate automatic exposure.

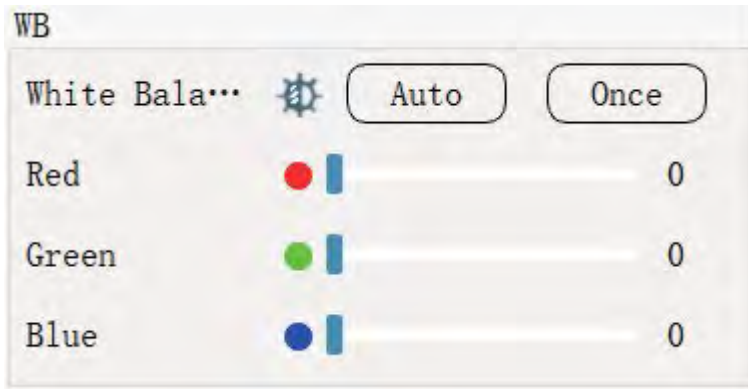


**Auto** automatic white balance is activated.

Manual adjustment of the red, green and blue levels is not possible 


The RGB values are regulated automatically by the camera.

**Auto** Automatic white balance is deactivated.


The parameters for red, green and blue levels can be adjusted manually. To do this, drag the slider to the left or  right accordingly.

**Once** This function is a one-time white balance. Click once and the camera automatically sets the values for red, green and blue levels.


### 10.3 Adjust parameters

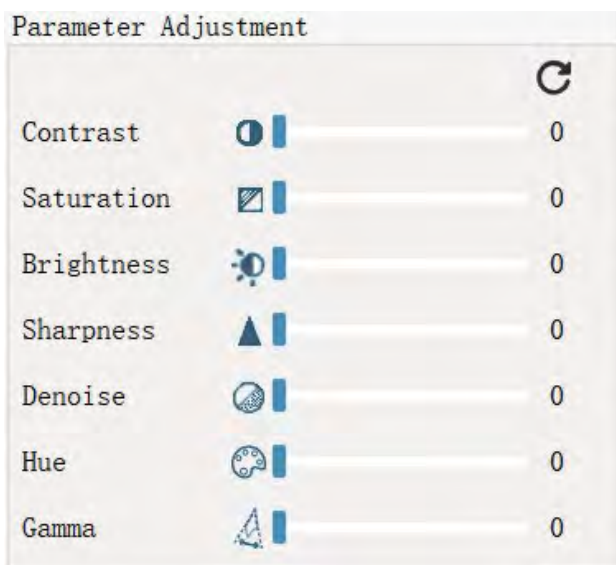


Drag the slider to the left or right to increase or decrease the values.


 Click on this button to reset all parameters to their default values.

## 10.4 Control elements

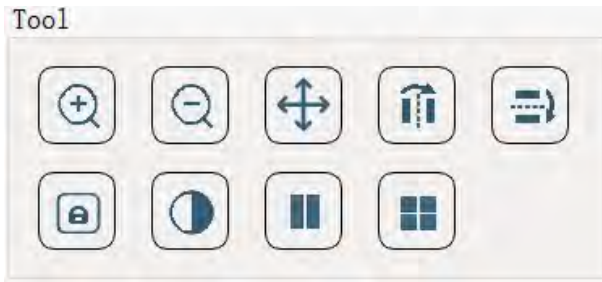



 **Magnification:** Magnification of the camera. The middle of the magnification is determined by the position of the previous magnification of the scroll wheel; if the scroll wheel has not been magnified, it is set to the middle magnification by default.

*Note:* The magnification is reset to 1.0X after opening the gallery. The magnification can also be increased using the mouse wheel.

 **Reduce:** Reduce the magnification of the camera. The center of the magnification is determined by the position of the last magnification, the next magnification is determined by the position of the last magnification, etc.

*Note:* The magnification is reset to 1.0X after opening the gallery. The mouse wheel can also reduce magnification.

 **Move:** If the current magnification factor is greater than 1.0, you can press and drag the left mouse button after clicking to move the current camera image.

 **Mirror horizontally:** The entire image is mirrored horizontally. The result of this process is that the point in the top right-hand corner swaps position with the point in the top left-hand corner and the point in the bottom right-hand corner swaps position with the point in the bottom left-hand corner.



Vertical mirroring: The entire image is mirrored vertically. The result of this process is that the dot in the top right-hand corner swaps position with the dot in the bottom right-hand corner and the dot in the top left-hand corner swaps position with the dot in the bottom left-hand corner.

To prevent damage to the focusing system, the left and right adjustment wheels of the coarse and fine drive must never be turned in opposite directions at the same time.



Freeze: Freezes the image displayed by the current camera.

*Note*: Increase in power, decrease in magnification, horizontal tilting, vertical tilting,

Monochrome and quad split screen cannot be performed in the current state. Freezing is automatically switched off when the gallery is opened.



Monochromatic: The color of the current frame is removed.



Comparison: After clicking, first open the gallery and then the photo or video to display the left and right frame. The left frame is the image of the current camera, and the right frame is the selected photo or video. The slider at the bottom moves to the left and right to show the incomplete display area of the camera image or the photo or video.


*Note*: The video contrast is only effective when the video is being played.



Four-part screen: After clicking, the screen is divided into four parts. Left-click on one of the screens to freeze it and click again to unfreeze it. Click on the quartered screen again to exit the quad screen. Exit the quad screen automatically when you open the gallery or the contrast.

*Note*: The measurement cannot be carried out if the quartering mode is open.


## 10.5 Measurements

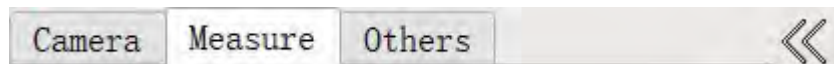


Click on **Measure** to switch to the measurement page, take measurements and change measurement attributes. As shown in the illustration on the right.

## 10.6 Measuring tools




 **Point:** Move the mouse pointer to the position to be measured, click the left mouse button and raise it to complete the drawing and display the corresponding coordinates.

*Note:* The coordinate values change depending on the position of the scale and the size of the calibration.

 **Line:** The first click selects the starting point, the second click completes the line drawing and displays the line and its length.

*Note:* The length of the line depends on the size of the calibration.

 **Parallel line:** The first click selects the starting point, the second click identifies a line, and the third click draws parallel and vertical lines between the two lines.

*Note:* The distance between the parallel lines depends on the size of the calibration.







Polyline: Click the left mouse button to get the different activities, after completion, click the right mouse button to stop, and show the total length of all lines.

*Note*: The length of the dashed line depends on the size of the calibration.



Arc: Select three points at different positions and draws an arc according to the distance between the three points and displays the radius and angle of the arc.

*Note*: The radius of the arc depends on the size of the calibration.



Radius circle: select the center of the circle, click the left mouse button and lift off, move the mouse, reach the point on the arc, click the left mouse button again, the operation is completed, select a circle corresponding to the two points, display the radius.

*Note*: The radius of the circle depends on the size of the calibration.



Diameter of the circle: Select a point on the arc, click the left mouse button, move to another point on the arc, click the left button again, the operation is completed, after the two points mark a circle, show the radius.

*Note*: The radius of the circle depends on the size of the calibration.



Three-point circle: With the left mouse button click three points in different positions, after the three points mark a circle, and show the radius of the circle.

*Note*: The radius of the circle depends on the size of the calibration.



Circle to line: the left mouse button click on two points in different positions, after the two points mark a circle, click the left button again, the point and the center of the circle mark a line after the point to do the vertical line of the line, show the distance from the point to the center of the circle.

*Note*: The distance between the circle and the line depends on the size of the calibration.



Concentric circles: Select the center of the circle, left-click and raise it, move the mouse to the point on the arc, click again, select the first circle corresponding to the two points, move the mouse to the end, left-click to raise it, select the second circle corresponding to the center of the circle and the end point, and display the radius of the double circle.

*Note*: The radius of the circle depends on the size of the calibration.



Double circle: Click three points at different positions with the left mouse button, select a circle corresponding to the three points, select three points at different positions again, click the left mouse button, select a circle in the same way, display the distance between the two circle centers.

*Note*: The length between the two centers depends on the size of the calibration.



Angle: The left mouse button click on three different positions of the point, after the three points to form two lines, show the angle of the two lines.



Vertical line: click the first time to select the starting point, click the second time to define a line, click the third time after the third point to form a vertical line, display a short line perpendicular to the line, display the length of the vertical line.

*Note*: The length of the vertical line depends on the size of the calibration.



Rectangle: Select the upper left corner of the rectangle, click the left mouse button to raise it, move the mouse to the lower right corner of the rectangle, click the left mouse button again to raise it, the operation is completed, the area and perimeter of the rectangle are displayed.

*Note*: The area and circumference of the rectangle change depending on the size of the calibration.



**Polygon**: Click the left mouse button to get each active point, after completion, click the right mouse button to stop, show the area and perimeter of the polygon.

*Note*: The area and circumference of the polygon change depending on the size of the calibration.



**Any line**: Move the mouse pointer to the position to be measured. Press the left mouse button and move it to draw any line. Press the left mouse button to complete the drawing.



**Text**: Move the mouse to the position to be measured, click on two points with the left mouse button, draw an arrow corresponding to the two points and open the input field. Click OK after the input in the input field is completed, and the entered text will be displayed from the end point of the input text. If you click Cancel, no text is entered, only the arrow is displayed.



**Move**: After clicking this icon, the mouse can capture the drawn points, line segments and circles. If the point clicked on is a measured point that has been drawn, you can change the size and position of the point. If this is not the case, switch to full movement.

*Note*: Polylines, polygons and arbitrary lines cannot move as a whole. If the whole thing moves, the text moves with it. If the position of the click is text, the text can be moved separately.



**Delete**: Click on this icon and left-click on the drawn measurement to delete the measurement.



**Delete all**: After clicking on this symbol, a reminder message appears. If yes, delete everything. If no, delete everything.



**Screenshot**: Capture the current image. Take screenshots without adding them to the toolbar with various prompts.



**Table:** To view the details of all currently drawn line segments, click Export to export the Excel file saved in the MEASUREMENT folder on the USB stick.

## 10.7 Calibration setting

Select the currently saved calibration options.



**Add calibration:** Click to enter the "Add Calibration" mode, enter the name and length, and then click the "Save" button to save the newly created calibration.



**Delete calibration:** Deletes the currently selected calibration.

*Note:* The first calibration cannot be deleted.



**Edit calibration:** Edit the currently selected calibration.



**Straight line calibration:** Switches the current calibration to a straight line.
















**Vertical calibration:** Switches the current calibration to vertical.







**Three-point circle calibration:** Switches the current calibration to three-point circle.

**Unit:** Switches the current calibration unit.

Add a new calibration: Click on  to add a new calibration, the delete calibration cannot be clicked at this point, the selection field for the calibration is empty. Select the calibration style    and click on the name input field, the keyboard appears. Enter the name or click  to close the keyboard and finish entering the name. Click on the input field for the length to call up the numeric keypad. After entering, click on  to close the numeric keypad and complete the name entry. Click on the drop-down field for the unit to complete the unit change. Click on  to save the calibration. Adding the calibration completed.

Edit the order of calibration: Click  select calibration, then the calibration selection bar, new calibration and delete calibration cannot be clicked. Select the desired calibration style, click on    the name and the keyboard will appear. After completing the entry, click ENTER or close  the keyboard and finish entering the name. Click on the input field for the length to call up the numeric keypad. When you have finished entering, click on to close  the numeric keypad and complete the name entry.

Click on the  drop-down field for the unit to change the unit and click on  to save the calibration. Edit calibration completed.

Delete calibration sequence: Select the calibration you want to delete in the selection bar  and click on  to delete the calibration.

## 10.8 Measurement attributes

Measurement Attribute			
Line Color	<input type="color" value="red"/>	Line Width	1
Font Color	<input type="color" value="red"/>	Font Size	10
Precision	0.1	Capture	N/A

Line Color (Line Color): After clicking, a drop-down field appears. After selecting the color, you can change the color.

*Note*: The measured color that was drawn is also changed.

Line width (Line Width): After clicking, a drop-down field appears. Select the width and change the line width.

*Note*: The width of the measured line width that was drawn is also changed.

Font color (Font Color): After clicking, a drop-down field appears. Select the color and change the font color.

*Note*: The font color of the drawn measurements is also changed.


Font size (Font Size): After clicking, a drop-down field appears in which you can select the size and change the font size.

*Note*: The font size of the drawn dimensions also changes.

Accuracy (Precision): The drop-down box appears after clicking, select the precision option and change the precision.

*Note*: The accuracy of the plotted measurements will also change.

## 10.9 Scale ruler

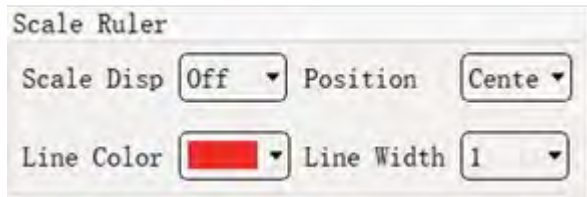


Display of the scale (Scale Disp): If you switch to On, the scale is displayed; if you switch to Off, the scale is closed.

Position: Switches the position displayed on the scale. If it is switched to the center, the reference point of the scale is in the center; if it is switched to the top left, the reference point of the scale is at the top left; if it is switched to the bottom right, the reference point of the scale is at the bottom right.

Color (Line Color): Click on the drop-down box that appears and select the color you want to change.

Width (Line Width): Switches the width of the scale of the scale.

## 10.10 Further settings

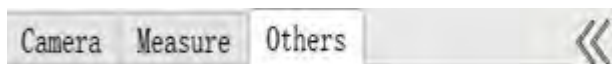

Click **Others to switch** to other pages to view and change the properties of the system.

## 10.11 Languages



## 10.12 Storage space

If the USB flash drive has several drives, you can select the drive on which the files are to be saved.

## 10.13 Resetting the data to the factory settings

After the factory settings have been restored, the exposure parameters are reset to the default values and automatic exposure is activated. Reset the white balance parameters to the default values and activate manual white balance. Contrast, saturation, brightness, sharpness, noise reduction, hue and gamma are reset to the default values.

Magnification zoom, screen movement, screen flip, freeze, monochrome, four-screen does not participate in factory restore.

The color of the measurement attribute is reset to red, the line width is reset to 2, the font color is reset to red, the font size is reset to 10, the accuracy is reset to 0.01 and the detection is reset to 10px.

The group color of the grid line V1-H1 is reset to red, the line is set to Off and the line width is reset to 1.

The group color of grid line V2-H2 is reset to orange, the line is reset to Off and the line width is reset to 1.

The group color of grid line V3-H3 is reset to yellow, the line is reset to Off and the line width is reset to 1.

The group color of the grid lines V4-H4 is reset to green, the lines are reset to Off and the line width is reset to 1.

The group color of grid line V5-H5 is reset to cyan, the line is set to Off and the line width is reset to 1.

The group color of the grid lines V6-H6 is reset to blue, the line is reset to Off and the line width is reset to 1.



The group color of the grid lines V7-H7 is reset to violet, the line is set to Off, the line width is reset to 1.

The color of the VB-HB group line is reset to white, the line to Off and the line width to 1.

The Group name option is reset to the first group V1-H1 and the display is reset to On. The scale display is reset to Off, the position is reset to the middle, the color is reset to red and the width is reset to 1.



## 10.14 Set date / time

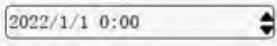
Change the time in  Click the mouse to change the year or month or day or hour or minute or second, the mouse wheel increases or decreases after the cursor flickers at the top, also you can click the small arrow on the right to change the setting system time. Click on  to modify the system time.






## 10.15 Switch between two displays


If you connect the monitor, you can control the display on the monitor or on the small screen here. If you connect HDMI before switching on, the system automatically displays the screen.

*Note:* You cannot switch when you are recording the screen. After you have clicked on the switch, the time-lapse photo will end automatically.

## 10.16 Move and hide toolbar

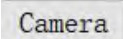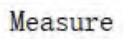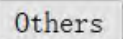 This button hides the side toolbar. If the current toolbar page is the measurement page.


If you click on the "Measurement tool"    button in the left toolbar (with the exception of "Delete all", "Screenshot" and "Table") and move the mouse over the left toolbar, the left toolbar is automatically hidden. If you want to show it, you can click the button  in the upper left corner to show it or click  to show the left status bar.

There is a move button  at the bottom of the camera, measurement and other pages in the left status bar, which can be clicked to move the left status bar to the right, and the simple measurement tool on the right will be moved to the left when it is displayed.


## 11 Toolbar




### 11.1 Camera

When the system recognizes that the USB flash drive is inserted and the available capacity is sufficient, click  to start recording. The top right corner will show that the photo is being saved. During this time, the camera button cannot be clicked again. The camera button can only be clicked when the message "Recording successful" is displayed in the top right-hand corner. The captured photo is named according to the current system time. Once the photos have been successfully saved, you can open the gallery to view the photos taken.

### 11.2 Camera timer

 If the system recognizes that the USB flash drive is inserted and the available capacity is sufficient, the timer page of the camera is displayed.

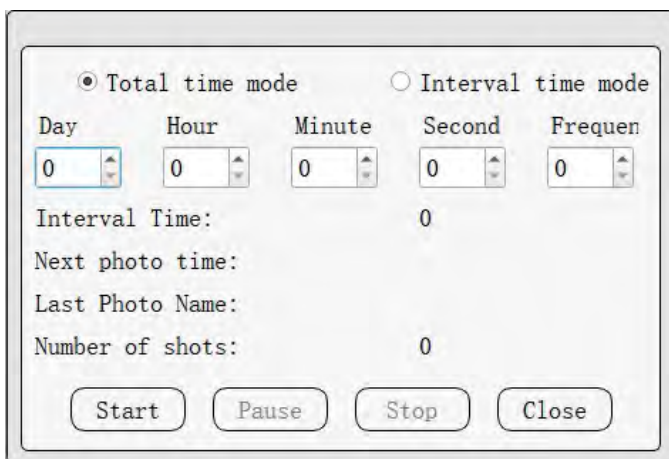A screenshot of a camera timer control window. It features two radio buttons at the top: "Total time mode" (selected) and "Interval time mode". Below these are five spinners for "Day", "Hour", "Minute", "Second", and "Frequen", each with a value of 0. Further down are labels for "Interval Time:" (0), "Next photo time:", "Last Photo Name:", and "Number of shots:" (0). At the bottom are four buttons: "Start", "Pause", "Stop", and "Close".

Total time mode (Total time mode): Set the timer mode to the number of photos taken in total time.

Interval time mode (interval tome mode): Set the interval at which a photo is taken.

Day (Day): Set the number of days. Click with the mouse, scroll wheel or click on the small button to control the increase or decrease.

Hour (Hour): Set the hours. Click with the mouse, scroll wheel or click on the small control button to increase or decrease.

Minute: Set the minutes. Click the mouse, scroll wheel or click the small button to increase or decrease the control.

Second (Second): Set the second speed. Click the mouse, scroll wheel or click the small control button to increase or decrease the speed.

Frequency: Specify the number of recordings. Click the mouse, scroll wheel or click the small control button to increase or decrease the frequency.

Interval Time (Interval Time): Displays the time interval for each photo taken in the current mode.

Time for next photo (Next photo time): Shows how much time is left until the next photo is taken.

Name of the last photo (Last photo name): Displays the name of the last photo taken.

Number of shots (Numbers of shots): Number of timed photos taken this time

Start: The camera timer starts.


*Note*: Time photography cannot be activated if the interval is less than two seconds.

Pause: The camera timer pauses.



Stop: Camera timer is stopped

Close: Closing the camera timer window and stopping the camera timer

### 11.3 Video recording

If the system recognizes that the USB flash drive is inserted and the available capacity is sufficient, click on the Video button  to start recording. The status of the button


is 

The recording time and the  indicator are displayed in the top right or left corner. Click the video button again to end the video recording. At this time, the upper right or left corner will display that the recording is being saved and the video recording button cannot be clicked. When the saving process is complete, the success of the recording and the. The recording button is  restored. You can resume recording at this point.

If the storage space on the USB flash drive is less than 50 MB, the recording will stop automatically. After saving, recording is automatically activated and continued until the storage space on the USB flash drive is less than 50 MB.

*Note:* If you disconnect the USB flash drive during recording, the video file may be damaged and can no longer be opened.


### 11.4 Gallery

If the system recognizes that the USB stick is inserted and is not recording, you can call up the gallery to view  photos and videos.



Switch to the image preview.



Switch to the video preview. A  appears above each video to remind you that the current file is a video file.



to the first page




to the last page




to the previous page




to the next page

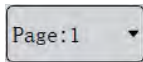

Displays the current page number. After clicking, a drop-down field appears from which you can select the page you want to jump to.

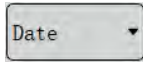

Sort order, date or name.

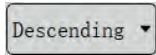

Sort order, ascending or descending.



Select an image or video and click Delete to delete the image or video. You can also select multiple items and then click Delete to delete them from the batch.




Home, Leave gallery.

## 11.5 Rename image and video files



After you have clicked on the image or video file, the selected image or video appears in the selected style.

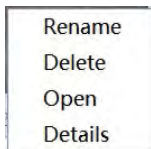

Right-click on the file to open the menu and select "Rename" or double-click on the file name to be renamed.



The cursor appears at the file name and calls up the keyboard.

After you have entered the name, click on





or close the keyboard




and complete the name change.


## 11.6 Delete image and video files


You can select one or more image/video files for deletion. After you have selected the file(s), the file(s) will have a blue  background. Then click on  to delete the file(s). Deleting the menu bar, which is called up by right-clicking on the file, can only delete the currently selected individual file.

## 11.7 Open image files








Double-click on the image to open it or right-click on the file to call up the menu and click on "Open" to open the image. After opening, the name of the current file is displayed in the control bar.

 to switch to the previous screen


 to switch to the next picture

 Leave gallery

## 11.8 Open video files

Double-click on the video to open it or right-click on the file, call up the menu and click on "Open" to open the video. After opening, the name of the current file is displayed in the control bar. After the video has been opened, it will play automatically, and the "Play" button will change to  Click  again to pause the current video and use the "Play" button to change to  Click  to switch to the previous video and play it automatically. Click on  to switch to the next video and play it automatically. Click to  pause the current video. Click  to exit the gallery.

## 11.9 File details

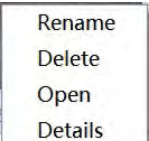 Right-click on the desired file, call up the menu bar and click on "Details" to display the details of the current file.

## 12 Optional accessories

<b>Bezeichnung</b>	<b>Artikel-Nr.</b>
Auxiliary objective 0.3x	OBB-A3240
Auxiliary objective 0.5x	OBB-A3225
Auxiliary objective 2.0x	OBB-A3226
Solder protection lens	OBB-A3227
Cleaning set for microscopes	OCS 901

## 13 Troubleshooting

Problem	Cause	Solution
Stains or dust on the picture	Impurities are on the sample	Clean sample
The image is blurred	Contamination is on the lens surface	Clean lens
	The focus is not correct	Adjusting the focus
The video image is not displayed	The external power cable is not connected	Connect the external power cable
The picture is too bright or too dark	The LED illumination is set incorrectly	Set the brightness of the LED ring light.
	Check the display settings	Open the "Settings" menu to set the brightness
The display does not respond	Overlaps in the system (too many operations) lead to error loops	Restart the system
The LED ring light does not light up when it is switched on	No power	Check the connection of the power cable
	The LED light has burnt out	Replace it with new lighting
The LED light has suddenly burnt out	The voltage is too high	Use a suitable external mains adapter
The illumination brightness is not sufficient	The voltage is too low	Use a suitable external mains adapter



## **14 Service**

If, despite studying these operating instructions, you still have questions about commissioning or operation, or if, contrary to expectations, a problem should occur, please contact your specialist dealer. The device may only be opened by trained persons and service technicians authorised by KERN.

## 15 Power supply

### 15.1 Opening video files



The microscope may only be connected to the mains supply if the specifications on the microscope (sticker) and the local mains voltage are identical.



**Important:**

- Check the mains cable for damage before commissioning
- Ensure that the power supply unit does not come into contact with liquids
- The mains plug must be accessible at all times.

## 16 Maintenance, servicing and disposal



Disconnect the appliance from the power supply before carrying out any maintenance, cleaning or repair work.

### 16.1 Maintenance and cleaning

The appliance must always be kept clean, and dust must be removed regularly. Clean the device with a damp, soft, lint-free cloth. Ensure that no moisture penetrates the housing. Do not use sprays, solvents, alcohol-based cleaners or abrasive cleaners, but only clear water to moisten the cloth.

Before wiping the appliance when wet, make sure that the power is switched off.

Glass components should preferably be wiped lightly with a lint-free cloth if they become dirty.

To wipe off oil stains or fingerprints from lens surfaces, moisten the lint-free cloth with a mixture of ether and alcohol (70/30 ratio) and then use it to clean the lenses

Ether and alcohol must always be handled with care, as they are highly flammable substances. It is therefore essential to keep them away from naked flames and electrical appliances that are switched on and off and only use them in well-ventilated rooms.

However, organic solutions of this type should not be used to clean other components of the appliance. This could cause changes to the paintwork. It is sufficient to use a neutral cleaning agent for this purpose.

Other cleaning agents for the optical components include

- Special cleaner for optical lenses
- Special optical cleaning cloths
- Bellows
- Brushes

If handled correctly and checked regularly, the microscope will function smoothly for many years.

However, should repairs be necessary, please contact your KERN dealer or our Technical Service.

### 16.2 Repair

Do not make any changes to the appliance and do not install any spare parts. Contact the manufacturer for repair or device inspection.

### 16.3 Disposal of the appliance and accessories



Old appliances and accessories should not be disposed of with household waste.

The operator must dispose of the packaging and appliance in accordance with the applicable national or regional legislation at the place of use.


The device consists of various components and materials, such as

- Electronic components (circuit boards, electrical cables)
- Plastic (housing)
- Metal (hook)

Improper disposal of the appliance can have harmful effects on people and the environment.

Proper and environmentally friendly disposal can prevent harmful effects and recover raw materials.

**Disposal of rechargeable batteries and batteries:**

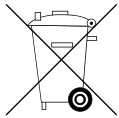

Rechargeable batteries and batteries do not belong in household waste.

The disposal of rechargeable batteries and batteries must be carried out by the operator in accordance with the applicable national or regional law of the place of use.

## 17 Battery Act

### Note according to the Battery Act - BattG:


#### INFORMATION




- The following information is valid for Germany.

In connection with the sale of batteries and rechargeable batteries, we are obliged as a dealer under the Battery Act to inform end users of the following:


- End users are legally obliged to return used batteries and rechargeable batteries.
- After use, batteries and rechargeable batteries can be returned free of charge to municipal collection centres or retailers. The normal end of use of the batteries/rechargeable batteries must be reached, otherwise precautions must be taken against short circuits.
- The return option is limited to batteries and rechargeable batteries of the type that we carry or have carried in our range and to the quantity that end users normally dispose of.
- A crossed-out wheelie bin means that you must not dispose of the batteries or rechargeable batteries in household waste under any circumstances. Old batteries or rechargeable batteries may contain harmful substances that can damage people and the environment if not disposed of correctly.

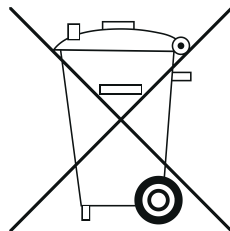


- Batteries containing harmful substances are labelled with a symbol consisting of a crossed-out dustbin and the chemical symbol (Cd = cadmium, Hg = mercury, or Pb = lead) of the heavy metal that is decisive for the classification as containing harmful substances.

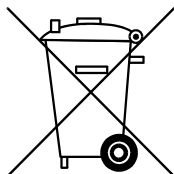

**Cd**

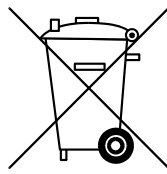

**Hg**

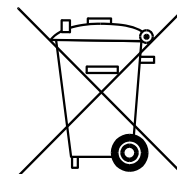

**Pb**

## 18 Further information

The illustrations may differ slightly from the product.

The descriptions and illustrations in these operating instructions are subject to change without notice. Further developments to the device may result in such changes.



All language versions include a non-binding translation.  
The original German document is binding.