

Operating Instructions | Mode d'emploi | Instrucciones de manejo

Original Operating Instructions | Mode d'emploi original | Traducción de las instrucciones de manejo originales

Cubis®

CUB Models: Microbalances and Ultramicrobalances

Modèles CUB: Microbalances et balances ultra-micro

Modelos CUB: Balanzas micro y balanzas ultramicro



SARTORIUS

imLab



www.imlab.eu - info@imlab.eu



+33(0)3 20 55 19 11



+32(0)16 73 55 72

Contents

1	About These Instructions	7
1.1	Validity.....	7
1.2	Target Groups.....	7
1.3	Symbols Used.....	7
1.3.1	Warnings in Operation Descriptions.....	7
1.3.2	Other Symbols Used.....	8
2	Safety Instructions	9
2.1	Intended Use.....	9
2.1.1	Modifications to the Device.....	9
2.1.2	Repairs to the Device.....	9
2.2	Qualifications of Personnel.....	10
2.3	Functionality of the Device Parts.....	10
2.4	Safety Information on the Device.....	10
2.5	Electrical Equipment.....	10
2.5.1	Damage to the Device's Electrical Equipment.....	10
2.5.2	Power Supply Unit and Power Supply Cable.....	10
2.6	Conduct in an Emergency.....	11
2.7	Accessories.....	11
2.8	Glass Breakage.....	11
3	Device Description	12
3.1	Device Overview.....	12
3.2	Draft Shield.....	13
3.3	Weighing Pan and Associated Components.....	14
3.4	Connections and Components on the Backside of the Weighing Module...	15
3.5	Connections on the Control Unit.....	15
3.6	Connections and Components on the Electronics Module.....	16
3.7	Weighing Compartment Lighting.....	17
3.8	Motorized Draft Shield.....	18
3.8.1	Sensor Area (not for Ultramicro- and Microbalances).....	18
3.9	Conformity-assessed Devices.....	18
3.10	Below-balance Weighing.....	18
4	Operating Design	19
4.1	Operating Elements in the Main Menu.....	19
4.2	Operating Elements in Task Management.....	20
4.3	Status Center.....	21
4.4	User Management.....	22
4.4.1	User Profiles.....	22
4.4.2	User Login.....	22
4.5	Weighing and Print Profiles.....	22
4.6	Applications and Tasks.....	22
4.7	Navigating the Menus.....	22

5	Installation	24
5.1	Scope of Delivery	24
5.2	Selecting an Installation Site	24
5.3	Unpacking	24
5.4	Preparing Below-balance Weighing	25
5.5	Installing a Microbalance with Glass Draft Shield	26
5.5.1	Positioning the Weighing Pan and Associated Components	26
5.6	Installing a Filter Balance with a Metal Ring Draft Shield	27
5.6.1	Positioning the Weighing Pan and Associated Components	27
5.6.2	Optimizing a Filter Weighing Pan with a Metal Ring Draft Shield for Left-handed Users	28
5.7	Acclimatizing	28
6	Getting Started	29
6.1	Connecting the Electronics Module	29
6.2	Connecting the Ethernet Cable	29
6.3	Installing the Power Supply Unit	30
6.4	Connecting the Power Supply	30
6.5	Connecting Accessories	31
6.6	Attaching Protective Caps and Covers	31
7	System Settings	32
7.1	Switching the Device On/Off and Activating Standby Mode	32
7.2	User Login or Logout	32
7.3	Performing System Settings	33
7.4	Using the Help Function	33
7.5	Activating Applications (QAPPs) and Adding to a Task	33
7.5.1	Activating Applications	33
7.5.2	Adding an Application to a Task	34
7.6	Switching Off the isoCAL Function	34
7.7	Configuring Motorized Opening and Closing of the Draft Shield (Only for Devices with a Motorized Draft Shield)	35
7.8	Adding Weighing and Print Profiles to a Task	35
7.9	Downloading Additional Information	35
7.10	Observing Warm-up Time	36
8	Operation	37
8.1	Opening and Closing the Motorized Draft Shield (Only for Devices with a Motorized Draft Shield)	37
8.1.1	Opening or Closing by Touching the Touch Sensors	37
8.2	Leveling the Device	37
8.3	Calibration, Adjustment or Linearization	37
8.3.1	Adjusting With the isoCAL Function	38
8.3.2	Internally Calibrating and Adjusting the Device	38
8.4	Preparing Weighings	39
8.5	Weighing	39
8.6	Running Application (Example)	40
8.6.1	Executing the "Unit change" Function	40

9	Cleaning and Maintenance	41
9.1	Preparing the Device for Cleaning.....	41
9.2	Cleaning the Device.....	41
9.3	Recommissioning.....	42
9.4	Performing a Software Update.....	42
9.5	Performing a QAPP Center Update.....	43
10	Malfunctions	44
10.1	Status Messages.....	44
10.2	Warning Messages.....	44
10.3	Troubleshooting.....	45
10.4	Operating Issues.....	45
11	Decommissioning	46
11.1	Decommissioning the Device.....	46
12	Transport	46
12.1	Transporting the Device.....	46
13	Storage and Shipping	47
13.1	Storing.....	47
13.2	Returning the Device and Parts.....	47
14	Disposal	48
14.1	Disposing of the Device and Parts.....	48
15	Technical Data	49
15.1	Dimensions and Weight.....	49
15.2	Ambient Conditions During Storage and Transport.....	49
15.3	Installation Conditions.....	49
15.3.1	Installation Site.....	49
15.3.2	Ambient Conditions at the Installation Site.....	50
15.4	Acclimatization Before Power Supply.....	50
15.5	Electrical Data.....	51
15.5.1	Power Supply.....	51
15.5.2	Safety of Electrical Equipment.....	51
15.5.3	Electromagnetic Compatibility.....	52
15.6	Warm-up Time to Reach Operating Temperature.....	52
15.7	Interfaces.....	52
15.7.1	Specifications for the COM-RS232 Interface.....	52
15.7.2	Specifications for the USB-A Interface.....	52
15.7.3	Specifications for the USB-B Interface.....	53
15.7.4	Specifications for the USB-C Interface.....	53
15.8	Recommended Calibration Weight.....	53
15.9	Conditions for isoCAL Function.....	53
15.10	Data Storage Device.....	54
15.11	Integrated Clock.....	54
15.12	Backup Battery.....	54
15.13	Materials.....	54

15.14	Cleaning Agents and Cleaning Procedures	55
15.14.1	Approved Cleaning Agents.....	55
15.14.2	Approved Cleaning Procedures	55
15.15	Metrological Data	56
15.15.1	Models CUB2.7S CUB10.6S CUB6.6S CUB3.6P.....	56
16	Accessories.....	57
16.1	Accessories	57
16.1.1	Printers and Communication.....	57
16.1.2	Hardware and Software for Pipette Calibration.....	58
16.1.3	Titanium Holder	58
16.1.4	Filter Balance and Antistatic Accessories.....	58
16.1.5	Special Applications	58
16.1.6	Weighing Tables.....	58
16.1.7	Weighing Accessories	59
17	Sartorius Service	60
18	Trademark Information.....	60
19	Conformity.....	60
19.1	EU Declaration of Conformity	60

1 About These Instructions

1.1 Validity

These instructions are part of the device; they must be read in full and kept in a safe place. These instructions apply to the device in the following versions:

Device	Model
Ultramicrobalance and microbalance	CUB10.6S-1XX-A CUB10.6S-1-XX-AC CUB6.6S-1XX-A CUB6.6S-1XX-AC CUB6.6S-1XX-F CUB6.6S-1XX-FC CUB3.6P 1XX-A CUB3.6P- 1XX-AC CUB2.7S 1XX-A CUB2.7S- 1XX-AC CUB2.7S 1XX-F CUB2.7S- 1XX-FC

1.2 Target Groups

These instructions are addressed to the following target groups. The target groups must possess the specified knowledge.

Target Group	Knowledge and Qualifications
Operator	The operator is familiar with the device and the associated work processes. The operator understands the hazards which may arise when working with the device and knows how to prevent them.*

* If a person in the target group operates the software interface of the device, they are also the "user".

1.3 Symbols Used

1.3.1 Warnings in Operation Descriptions

WARNING

Denotes a hazard that may result in death or serious injury if it is **not** avoided.

CAUTION

Denotes a hazard that may result in moderate or minor injury if it is **not** avoided.

NOTICE

Denotes a hazard that may result in property damage if it is **not** avoided.

1.3.2 Other Symbols Used

- ▶ Required action: Describes activities that must be carried out. The activities in the sequence must be carried out in succession.
- ▷ Result: Describes the result of the activities carried out.
- [] Refers to operating and display elements. Indicates status, warning, and error messages.
- M Indicates information for legal metrology for conformity-assessed (verified) devices. Conformity-assessed devices are also referred to as “verified” in these instructions.

Figures in these instructions

Depending on the device configuration, the figures depicting the device and operating display may differ slightly from the supplied device. The variants shown in these instructions are examples.

2 Safety Instructions

2.1 Intended Use

The device is a high-resolution balance, which can be used in laboratories. The device is intended to accurately determine the mass of materials in liquid, paste, powder, or solid form.

Appropriate vessels must be used for loading certain materials.

The device can be operated as follows:

- In individual operation
- Connected to a PC
- Integrated into a network

The device is intended solely for use in accordance with these instructions. Any other use is considered **improper** and can impair the protection supported by the device, e.g., protection against mechanical dangers.

Foreseeable misuse

The following applications are **not** permitted: Operation under any atmosphere other than the normal atmosphere.

Operating conditions for the device

Do **not** use the device in potentially explosive environments. Only use the device indoors.

Do **not** use physical measures to alter the as-delivered state of the device, and only connect approved accessories (see Chapter "16 Accessories", page 57).

The device may only be used with the equipment and under the operating conditions described in the Technical Data section of these instructions.

2.1.1 Modifications to the Device

If the device is modified: Persons may be put at risk. Device-specific documents and product approvals may lose their validity.

For queries regarding modifications to the device, contact Sartorius.

2.1.2 Repairs to the Device

Device repairs may only be carried out by persons with specialized knowledge of the device. If the device is **not** repaired by a specialist: Persons may be put at risk. Device-specific documents and product approvals may lose their validity.

Sartorius recommends that any repair work, even that carried out after the end of the warranty period, is carried out by Sartorius Service or after consulting with Sartorius Service.

2.2 Qualifications of Personnel

Persons who do not possess adequate knowledge about how to use the device may injure themselves and other persons.

If a particular qualification is required for an activity: The target group is specified. If **no** qualification is specified: The activity can be performed by the "Operator" target group.

2.3 Functionality of the Device Parts

Non-functioning device parts, e.g., as a result of damage or wear, can cause malfunctions. There is a risk of injury to persons.

- ▶ If device parts are **not** functioning: Do **not** use the device.

2.4 Safety Information on the Device

Symbols, e.g., warnings and safety stickers, are safety information for handling the device. Missing or illegible safety information may result in this information **not** being observed. There is a risk of injury to persons.

- ▶ Do **not** conceal, remove, or modify the symbols.
- ▶ Replace the symbols if they become illegible.

2.5 Electrical Equipment

2.5.1 Damage to the Device's Electrical Equipment

Damage to the device's electrical equipment, e.g., damaged insulation, can be life-threatening. Contact with live parts represents a danger to life.

- ▶ If the electrical equipment of the device is defective, disconnect the power supply and contact Sartorius Service.
- ▶ Keep live parts away from moisture. Moisture can cause short circuits.

2.5.2 Power Supply Unit and Power Supply Cable

The use of an **unauthorized** power supply unit or power supply cable may cause life-threatening injuries as a result of electric shocks, for example.

- ▶ Only use the original power supply unit and power supply cable.
- ▶ If the power supply unit or power supply cable needs to be replaced: Contact Sartorius Service. Do **not** repair or modify the power supply unit or power supply cable.

2.6 Conduct in an Emergency

If an emergency occurs, e.g., due to malfunctions of the device or dangerous situations: Persons may be injured. The device must be immediately taken out of operation:

- ▶ Disconnect the device from the power supply.
- ▶ Prevent the device from recommissioning.

2.7 Accessories

The use of unsuitable accessories can affect the functionality and operating reliability of the device and has the following consequences:

- Risk of injury to persons
 - Damage, malfunctions, or failure of the device
- ▶ Only use accessories that have been approved by Sartorius for this device.

2.8 Glass Breakage

Glass components can break if they fall or are handled incorrectly. Glass fragments can cause cuts.

- ▶ Do **not** use sharp or hard objects on the operating display.
- ▶ Do **not** allow items to fall onto the control unit.
- ▶ In the event of damage to the control unit or draft shield, do **not** use the device. Contact Sartorius Service.

3 Device Description

3.1 Device Overview



Fig. 1: Microbalance with motorized glass draft shield and electronics module (example)

Pos.	Name	Description
1	Weighing compartment	
2	Leveling foot	For manually leveling the device
3	Electronics module	
4	Manufacturer's ID label	<ul style="list-style-type: none"> – Contains additional information about the device, e.g., serial number and the metrological and technical data. – Not visible.
5	Operating display	Touchscreen
6	Touch sensor	<ul style="list-style-type: none"> – Only for microbalances with motorized draft shield: Opens and closes the draft shield. – Reacts to touch on the left and right on the top, and on the side panels. – If active: The symbol illuminates at half brightness. – If touched: The symbol illuminates at full brightness.
7	Control unit	
8	Weighing module	

3.2 Draft Shield



Fig. 2: Microbalance with motorized glass draft shield and filter microbalance with manual metal ring draft shield (example)

Pos.	Name	Description
1	Metal draft shield cover	Made from metal, with handle, can be removed.
2	Metal ring draft shield	Consists of 2 metal rings that have been placed into each other with a side opening, can be rotated manually.
3	Glass draft shield	Made from glass, with side opening, can be rotated manually or via a motor. Can be removed.
4	Sensor window	For the motion sensor, controls opening and closing of the motorized glass draft shield

3.3 Weighing Pan and Associated Components

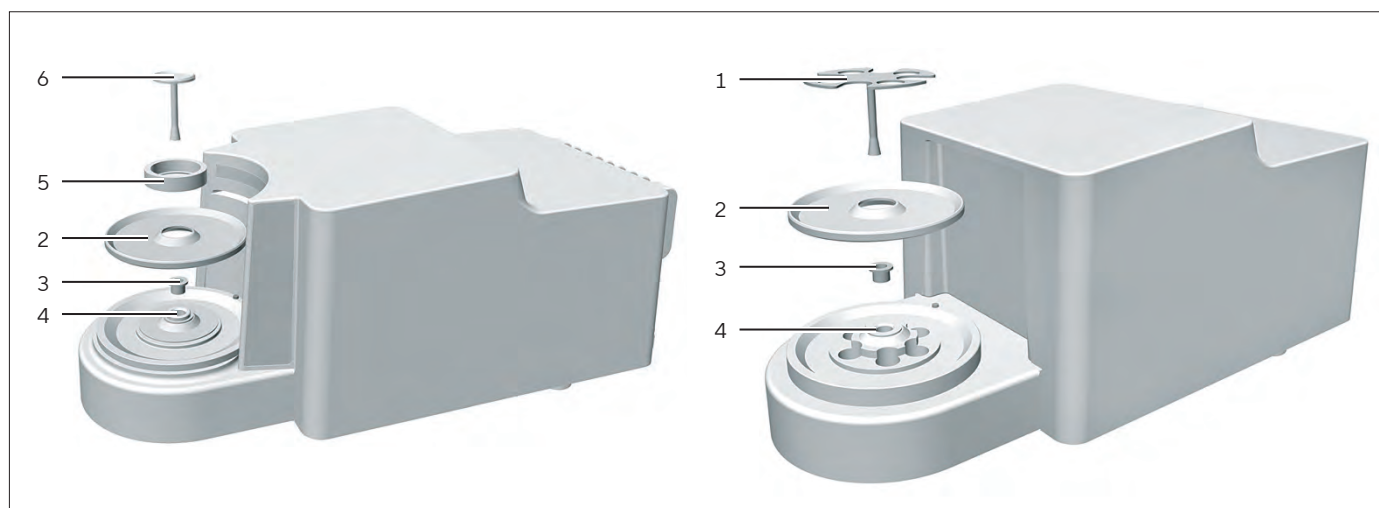


Fig. 3: Microbalance with motorized glass draft shield and filter microbalance with manual metal ring draft shield (example)

Pos.	Name	Description
1	Filter weighing pan	
2	Shield plate	
3	Connector	For models CUB2.7S... only
4	Pan retainer	
5	Internal draft shield	For models CUB2.7S...A only
6	Weighing pan	

3.4 Connections and Components on the Backside of the Weighing Module

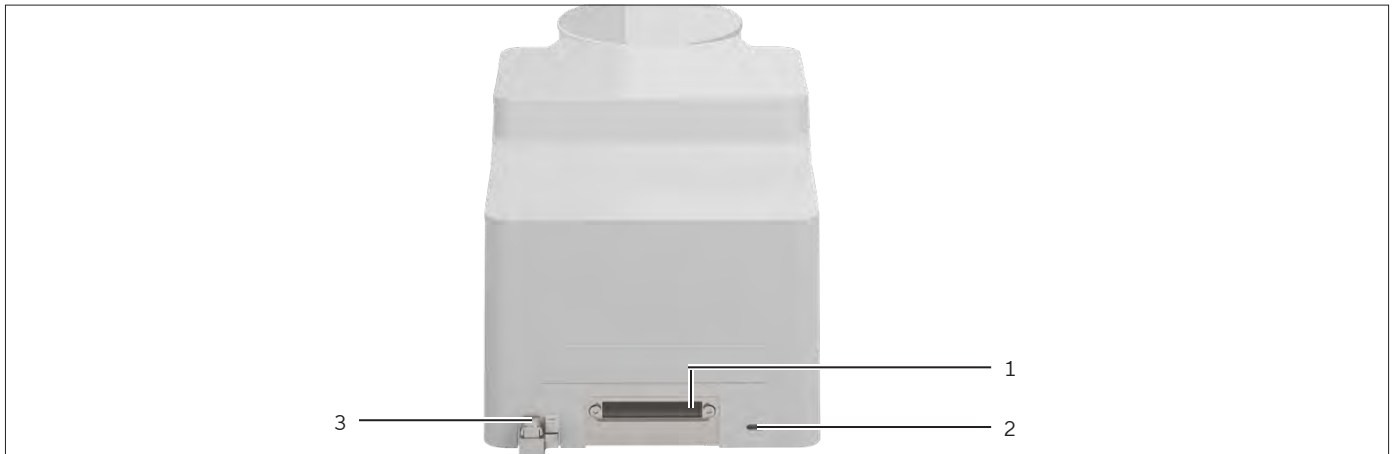


Fig. 4: Connections on the weighing module of the microbalance (example)

Pos.	Name	Description
1	Electronics module connection	For connecting the electronics module to the weighing module
2	Slot	For connecting a "Kensington" anti-theft device
3	Terminal	For connection of a potential equalization

3.5 Connections on the Control Unit

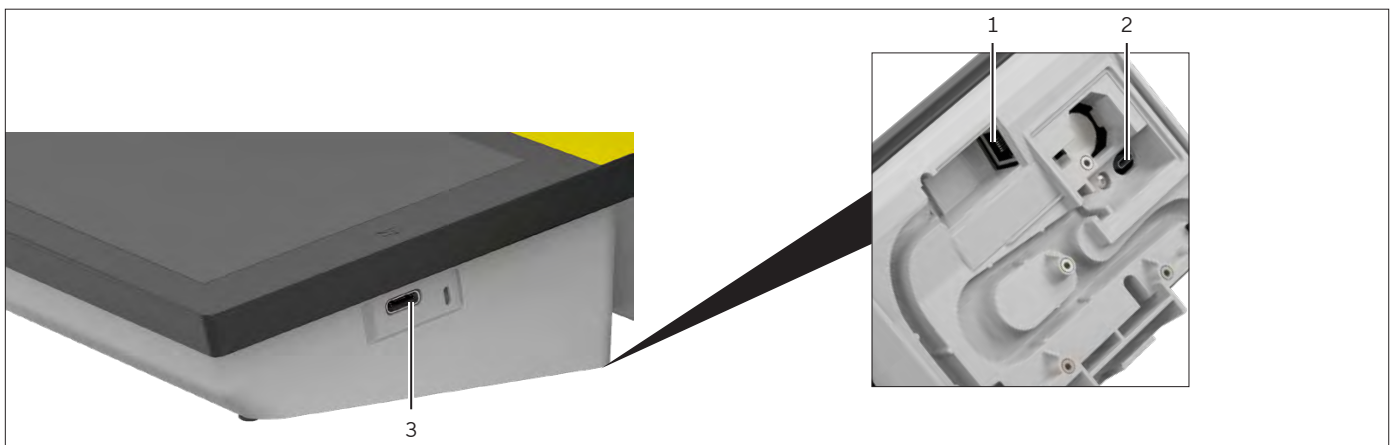


Fig. 5: Connections on the control unit (example)

Pos.	Name	Description
1	Weighing module connection	For connecting the device.
2	Ethernet connection	For connection to a network.
3	USB-C connection	For connecting accessories.

3.6 Connections and Components on the Electronics Module

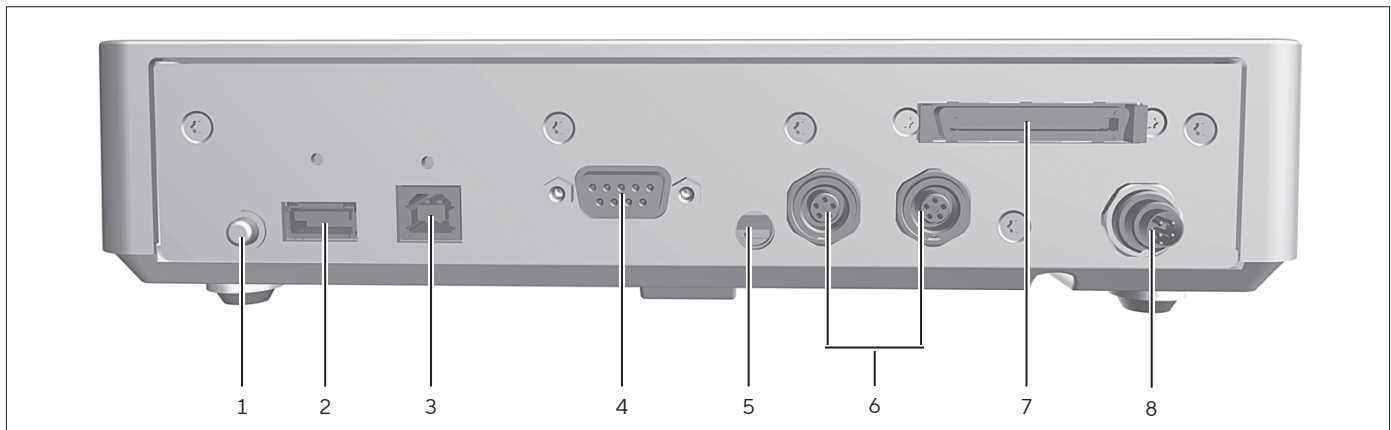


Fig. 6: Connections on the electronics module (example)

Pos.	Name	Description
1	On key	Switches the device back on again after switching off via the software.
2	USB-A connection	<ul style="list-style-type: none"> – For USB accessories, e.g., printers, USB mass storage devices, bar-code scanners. – Sealed with a plastic protective cap. – The protective cap is secured to the device.
3	USB-B connection	<ul style="list-style-type: none"> – For connection to a PC. – Sealed with a plastic protective cap. – The protective cap is secured to the device.
4	COM-RS232 connection	<ul style="list-style-type: none"> – 9-pin, for connection to a PC or PLC – Sealed with a plastic protective cap. – The protective cap is removable.
5	Access switch	<ul style="list-style-type: none"> – Protects the device from changes to the device settings. Is sealed for conformity-assessed devices. – Sealed with a plastic protective cap. – The protective cap is removable.
6	Peripheral connection	<ul style="list-style-type: none"> – For connecting Sartorius accessories. – Sealed with a plastic protective cap. – The protective cap is removable.
7	Weighing module connection	For connecting the electronics module to the weighing module.
8	Power supply	For connection to the power supply.

3.7 Weighing Compartment Lighting



Fig. 7: Weighing compartment lighting (example)

Pos.	Name	Description
1	Weighing compartment lighting	Intended for illuminating the weighing compartment. The intensity of the lighting can be adjusted in the menu. When the control unit displays an error message: The weighing compartment is illuminated in orange.

3.8 Motorized Draft Shield

3.8.1 Sensor Area (not for Ultramicro- and Microbalances)

The device has a motion sensor that opens the draft shield automatically as soon as a motion is detected within the sensor area. Following that the draft shield will be automatically closed again. The operation is optimized for left- and right-handed users due to the motion sensor.

The motion sensor has a right and a left sensor window.

The motion sensor has two modes:

- Proximity sensor (one-hand operation)
- Crossover method (two-hand operation)

The motion sensor sensitivity can be adjusted in steps.

3.9 Conformity-assessed Devices

Some settings of conformity-assessed models are protected against operator changes, e.g., "external adjustment". This measure is intended to ensure the suitability of the devices for use in legal metrology.

3.10 Below-balance Weighing

The device is suitable for below-balance weighing. Samples can be suspended for weighing using below-balance weighing, e.g., samples, which do **not** fit on the weighing pan. Below-balance weighing is possible under the following conditions:

- The device must be set up on a weighing table with recess.
- For below-balance weighing, a below-balance weighing hook must be inserted in the device base. The below-balance weighing hook is available as an accessory (see Chapter "16 Accessories", page 57).

M

In legal metrology:

- The below-balance weighing equipment may **not** be used.
- The cover of the below-balance weighing may **not** be removed.

4 Operating Design

4.1 Operating Elements in the Main Menu

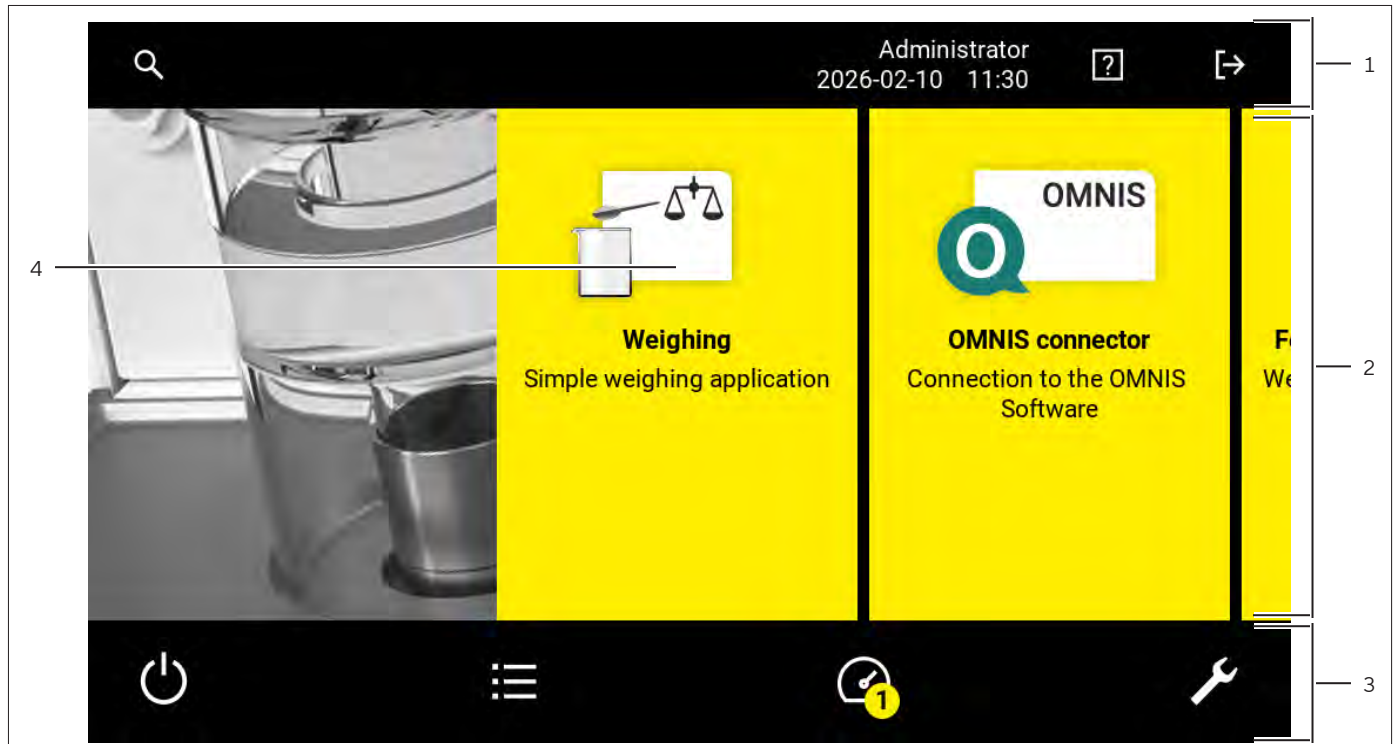


Fig. 8: Operating elements in the main menu (example)

Pos.	Name	Description
1.	Navigation and function bar	<ul style="list-style-type: none"> – Enables navigation and searching in menus and lists. – In the “Settings” menu: Displays the name of the menu.
2.	Available tasks	Displays all tasks available for the active user.
3.	Function bar	Displays available submenus and operating functions for the current display and current user.
4.	Task	Starts the described task.

4.2 Operating Elements in Task Management

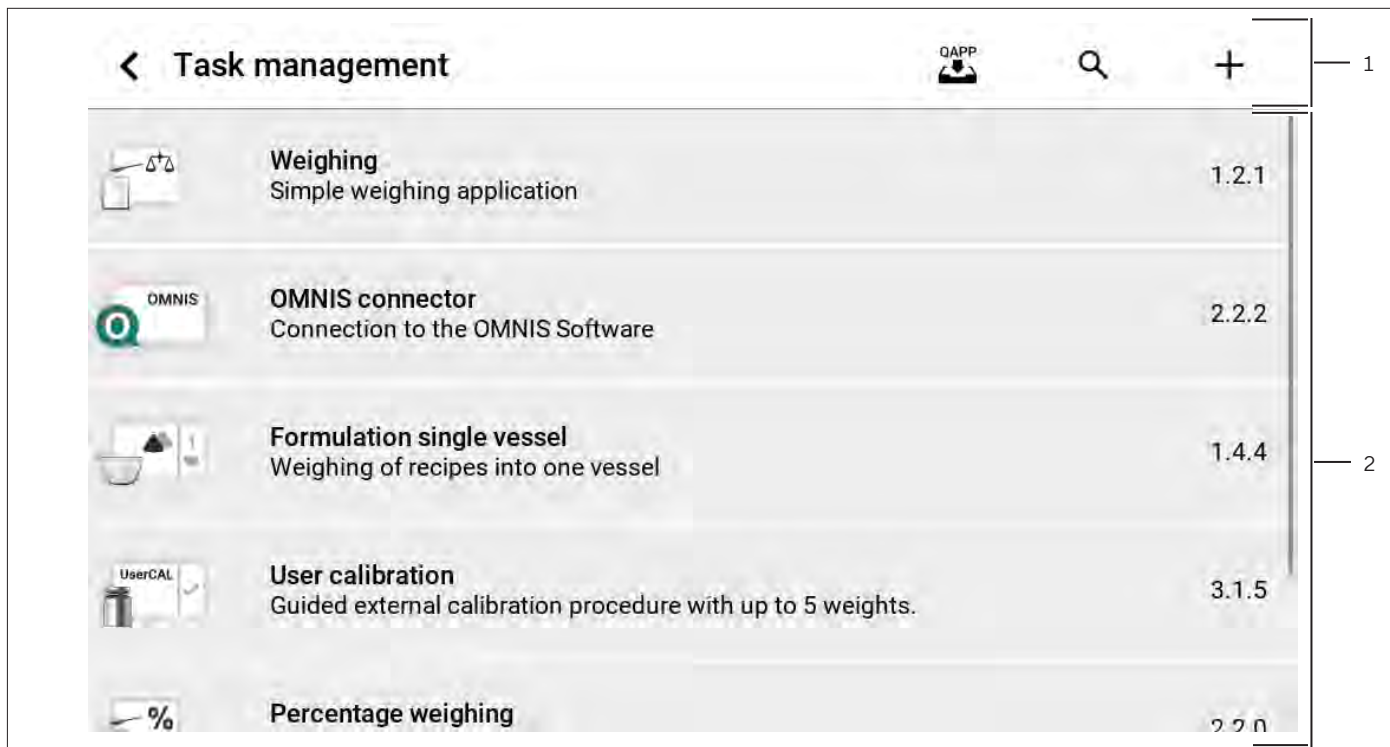


Fig. 9: Operating elements in task management (example)

Pos.	Name	Description
1.	Navigation and function bar	<ul style="list-style-type: none"> – Enables navigation and searching in menus and lists. – Enables the addition of tasks. – Opens the QAPP Center. – Displays the name of the menu.
2.	Available tasks	<ul style="list-style-type: none"> – Displays all available tasks. – Opens a summary of the properties for the displayed task.

4.3 Status Center

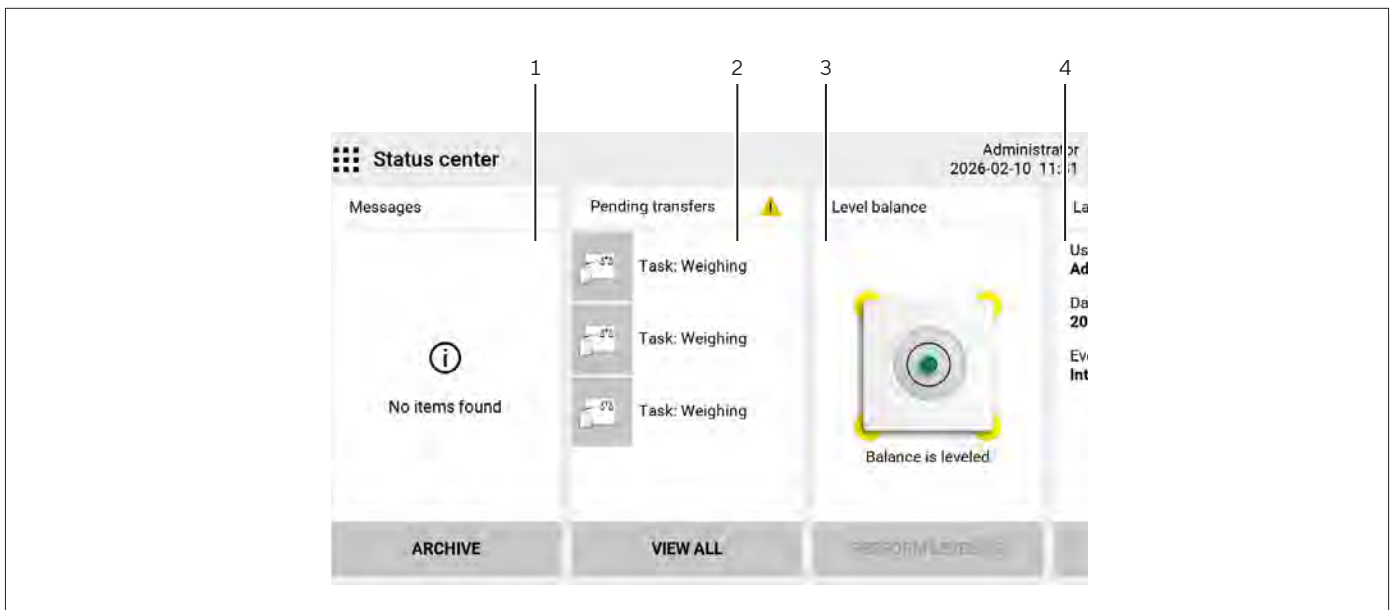


Fig. 10: Status center (example)

Pos.	Name	Description
1	Messages	Displays information, warning, and error messages.
2.	Pending transfer	If there are pending transfers: Displays the pending transfers for completed tasks.
3.	Leveling status	Displays the status of the level, starts leveling.
4.	Further categories	Displays further categories that are made visible when swiping to the left, e.g., <ul style="list-style-type: none"> – Status for the device – Calibration and adjustment report – Audit trail

4.4 User Management

4.4.1 User Profiles

In the factory, four user profiles are created for the device. One role is assigned to each user profile. Each role has rights to operate the device. The rights assigned to each role depend on which device functions the user has to use. User profiles can be adapted.

Additional user profiles, roles, and rights may be activated in the device settings for devices with licensed user management.

4.4.2 User Login

The user must log in to the login display with a user profile. Various setting options and tasks are displayed in the operating display depending on the user profile and role.

4.5 Weighing and Print Profiles

Weighing and print profiles can be created. These profiles can be assigned to a task.

Preset profiles can be used in a task. The preset profiles can be adjusted individually to the application and saved in new profiles for weighing and printing.

4.6 Applications and Tasks

Applications (QAPP applications) are grouped together in QAPP packages. Depending on the model, the device can be supplied with some freely accessible applications. These applications can be used to carry out the most important functions.

Subsequent licensing (package QP99) covers all other applications and can be activated in the QAPP Center for a fee.


An application must be configured as a task before it can be used. Specific settings need to be applied for this, with the help of the assistant function. Tasks are visible to all users who have the required role to perform the task.

4.7 Navigating the Menus

Procedure

- ▶ To open an application or a menu from the main menu:
- ▶ Press the application button or required menu button in the function bar.
- ▷ The application or menu is opened and the name of the open menu is opened in the navigation bar.
- ▶ To return to the main menu from other displays: Press the [Back] or [Menu] button.




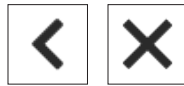

- ▶ To return to the next higher menu level: Press the [Back] or [Cancel] button.



- ▶ To scroll through the available tasks or categories in a horizontal menu, e.g., main menu or status center: Swipe to the left or right on the operating display.

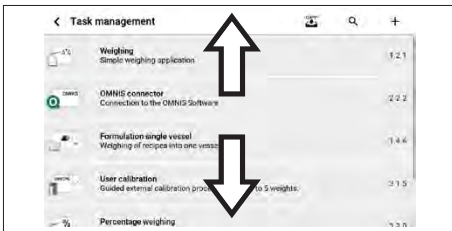


- ▶ To scroll through the lists in a vertical menu: Swipe the list downwards or upwards.




- ▶ If a value needs to be selected from a list:
 - ▶ Press the preferred value.
 - ▶ Confirm the selection with [OK].
- ▷ The selected value is saved and the list closes.

- ▶ If elements from a display need to be filtered or a display needs to be browsed:
 - ▶ Press the [Search] or [Filter] button.

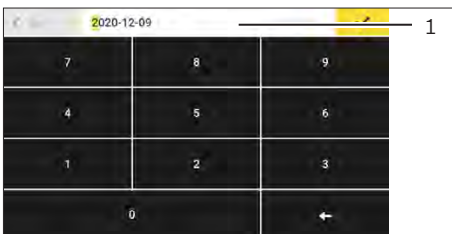

- ▷ The keypad is displayed.
 - ▶ Type the search value or value to be filtered into the input field (1) using the keypad and confirm with [OK].
- ▶ To close the input field without searching or filtering: Leave the input field empty and press the [OK] button.

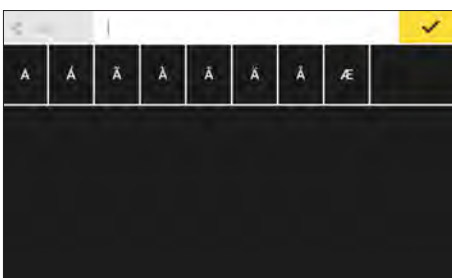

- ▶ If language-specific characters need to be entered using the keypad:
 - ▶ Press and hold a letter on the keypad.
 - ▷ If language-specific characters are available for the letter being pressed: The display containing the language-specific characters opens.
 - ▶ Press the required language-specific character.

5 Installation

5.1 Scope of Delivery

Item	Quantity
Device	1
Weighing pan	1
For filter microbalances, ultramicrobalances and microbalances: Filter weighing pan	1
For models CUB2.7S...: Internal draft shield	1
Shield plate	1
For models CUB2.7S...: Connector	1
Power supply unit	1
Country-specific power supply cable	1 – 4*
In-use dust cover for the control unit	1
Electronics module	1
Connection cable for electronics module	1
Operating Instructions	1 – 2*

* Quantity is country-specific.

5.2 Selecting an Installation Site

Procedure

- ▶ Ensure that the following conditions are met at the installation site (see Chapter “15.3 Installation Conditions”, page 49).
- ▶ **NOTICE** Risk of damage to the power supply unit from argon! Observe the instructions for using argon (see Chapter “15.3.2 Ambient Conditions at the Installation Site”, page 50).

5.3 Unpacking

Procedure

- ▶ Lift the device with the polystyrene out of the packaging.
- ▶ Place the device in the polystyrene on its side.
- ▶ Lift the polystyrene off the device.
- ▶ **⚠ CAUTION** Glass breakage due to the incorrect handling of the device! Do **not** lift the device by the draft shield. Only lift the device by its base.
- ▶ Place the device on its base.
- ▶ Keep all parts of the original packaging, e.g., to return the device.

5.4 Preparing Below-balance Weighing

The device can be configured for below-balance weighing. Samples can be suspended for weighing using below-balance weighing, e.g., samples, which do **not** fit on the weighing pan.

For below-balance weighing, the cover plate must be removed and the device set up on a weighing table with recess.

M

In legal metrology:



- The below-balance weighing equipment may **not** be used.
- The cover of the below-balance weighing equipment may **not** be opened.

Material: 1 draft protection shield
1 weighing table with recess

Tool: 1 Torx Allen key, T20

Procedure

- ▶ If a draft shield is mounted: Remove the glass draft shield.
- ▶ Remove the base plate of the weighing compartment.
- ▶ Place the device on its side.

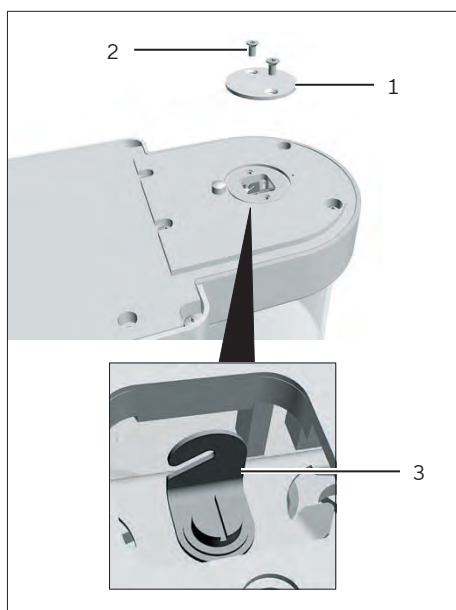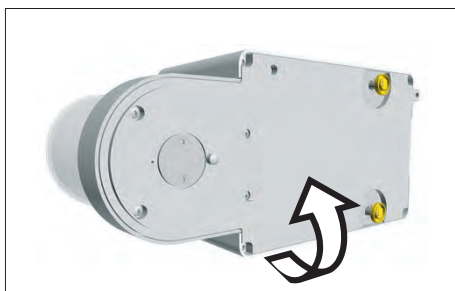


- ▶ Unscrew both screws (2) from beneath the weighing module with the Torx Allen key.
- ▶ Lift out the cover plate (1).
- ▶ Install a draft protection shield.

- ▶ Set up the device on the weighing table with recess. The hook for below-balance weighing (3) may **not** touch the weighing table.
- ▶ Install the draft protection shield.


5.5 Installing a Microbalance with Glass Draft Shield

5.5.1 Positioning the Weighing Pan and Associated Components


Procedure

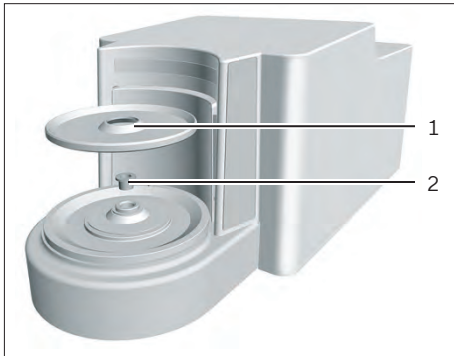

- ▶ For a device with a connector: Place the connector (2) on the base of the weighing compartment.
- ▶ Place the shield plate (1) on the base of the weighing compartment.

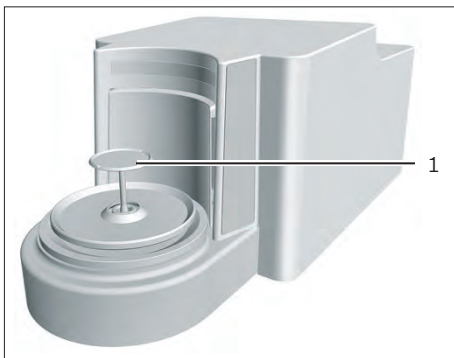

- ▶ Place the weighing pan (1) through the opening in the base plate of the weighing compartment into the pan retainer.
- ▶ To mount the weighing pan: Turn the weighing pan slightly while pressing down lightly.

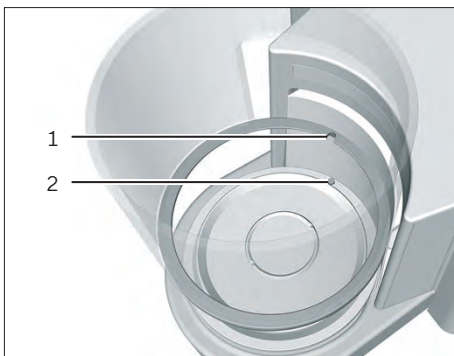



- ▶ For devices with an internal draft shield: Place the internal draft shield on a shield plate.
- ▶ Center and insert glass draft shield with recess (1) above the snap-in peg (2).

5.6 Installing a Filter Balance with a Metal Ring Draft Shield


5.6.1 Positioning the Weighing Pan and Associated Components

Procedure

- ▶ For a device with a connector: Place the connector on the base of the weighing compartment.
- ▶ Place the shield plate (2) on the base of the weighing compartment.
- ▶ Place the metal ring draft shield (1) on the base of the weighing compartment. The inner metal ring of the draft shield must be inserted into the outer metal ring.

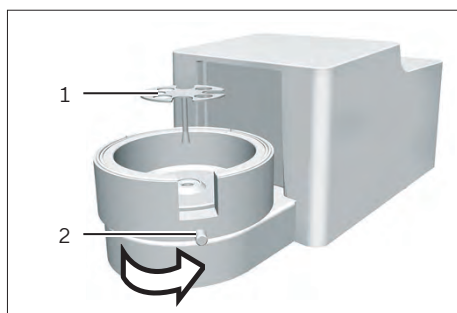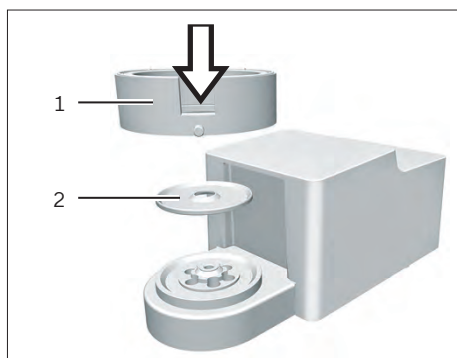

- ▶ Turn the opening of the metal ring draft shield in the direction of the arrows.
- ▶ To mount the metal ring draft shield: Turn the screw (2).
- ▶ Place the filter weighing pan (1) or the weighing pan through the opening in the shield plate into the pan retainer.
- ▶ **NOTICE** The filter weighing pan or weighing pan can be damaged due to improper assembly. To mount the inserted filter weighing pan or weighing pan: Screw in the filter weighing pan or weighing pan with slight pressure to correctly position it in the pan retainer.

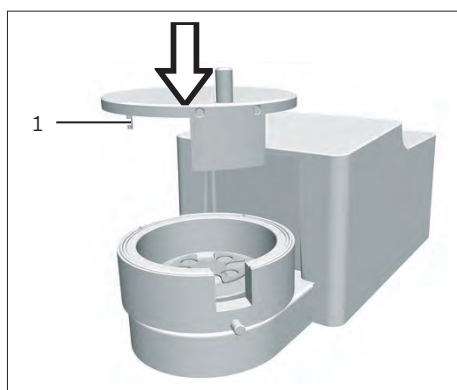

- ▶ Insert the bearing pin (1) in the draft shield cover into the recess of the draft shield.
- ▶ Place the draft shield cover on the draft shield.


5.6.2 Optimizing a Filter Weighing Pan with a Metal Ring Draft Shield for Left-handed Users

Requirements


The metal ring draft shield has been mounted.

Procedure

- ▶ Lift the draft shield cover from the metal ring draft shield and turn.
- ▶ Pull out the bearing pin (1) from the underside of the draft shield cover and insert into the opening on the opposite side.

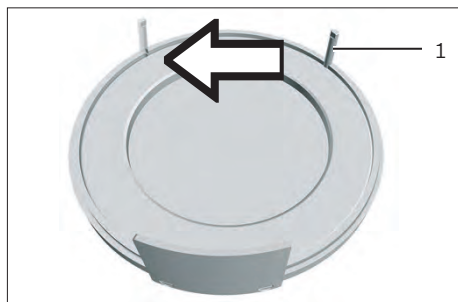

- ▶ Place the draft shield cover on the metal ring draft shield.
- ▶ Loosen the screw (1) on the metal ring draft shield and turn all parts of the draft shield in the direction of the arrows by 90°.
- ▶ Retighten the screw.

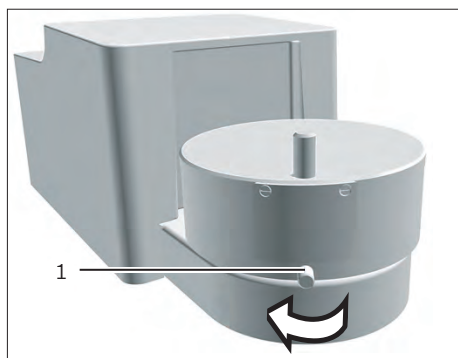

5.7 Acclimatizing

When a cold device is brought into a warm environment: The temperature difference can lead to condensation from humidity in the device (moisture formation). Moisture in the device can lead to malfunctions.

- ▶ Allow the device to acclimatize at the installation site (acclimatization duration see Chapter "15.4 Acclimatization Before Power Supply", page 50). The device must be disconnected from the power supply during that time.

6 Getting Started

Procedure

- ▶ **NOTICE** Improper connection may damage the device! If the device is connected to electronic components, e.g., printer, PC: The device must be disconnected from the power supply. Ensure that the device is disconnected from the power supply.
- ▶ Connect the device with the electronic components (see electronic components instructions).

6.1 Connecting the Electronics Module

Procedure




- ▶ Connect the connection cable to the electronics module's weighing module connection.
- ▶ Connect the other end of the connection cable to the weighing module's electronics module connection.
- ▶ To lock the connection cable: Lock the plugs of the connection cable onto both connections with two clicks in each case.
- ▶ **NOTICE** Damage to the device due to incorrect connection!
 - ▶ Check the correct fit of the plug contacts.
 - ▶ There should be **no** tension on the connection cable, e.g., do **not** install directly against a wall.

6.2 Connecting the Ethernet Cable

Material: 1 Ethernet cable


Procedure

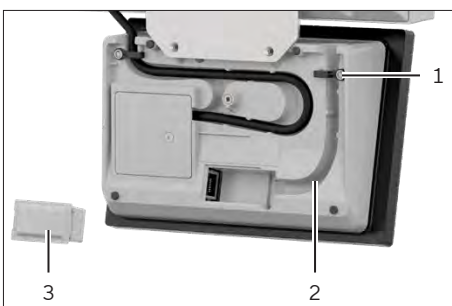


- ▶ Turn over the electronics module with the control unit.
 - ▶ Remove the cover (3) of the Ethernet connection on the underside of the control unit.
 - ▶ Plug the Ethernet cable into the Ethernet connection.
 - ▶ Place the Ethernet cable into the cable channel (2) and turn the cable lock (1) using the cable.
- ▶ Place the electronics module with control unit back on the device base on a level surface.

6.3 Installing the Power Supply Unit

Procedure

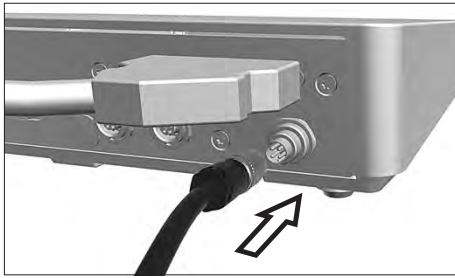

- ▶ Plug the DC supply cable of the AC adapter into the electronics module's "Power supply" connection and tighten the threaded fitting.

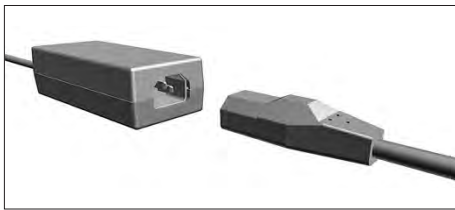

- ▶ Plug the power supply cable into the power supply unit connection.

6.4 Connecting the Power Supply

Requirements

The acclimatization time has been observed and the device has adjusted to room temperature (see Chapter 15.4, page 50).

Procedure

- ▶ Check whether the country-specific power plug matches the power supplies at the installation site.
 - ▶ If required: Contact Sartorius Service.
- ▶ **NOTICE** Damage to the device due to excessive input voltage! Check whether the voltage specifications on the power supply unit match those of the power supply at the installation site.
 - ▶ If the input voltage is too high or too low: Do **not** connect the device to the power supply.
 - ▶ Contact Sartorius Service.
- ▶ Connect the device to the power supply at the installation site. To do this, connect the power plug from the power supply cable into the power socket.
- ▶ The device is switched on and performs the initial functions for device startup.

6.5 Connecting Accessories

Accessories can be connected to the device.

Requirements

The accessories must be suitable for the device (see instructions for the accessories).

Procedure


- ▶ Ensure that the device is disconnected from the power supply.
 - ▶ If required: Disconnect the device from the power supply.
- ▶ Connect the accessories to the appropriate connections on the device (connections for accessory parts see instructions for accessories).
- ▶ Connect the device to the power supply.

6.6 Attaching Protective Caps and Covers

If connections of the device are **not** being used during operation: We recommend sealing the connections with the protective caps provided.

Procedure

- ▶ Check whether all unused connections have been sealed with a protective cap.
 - ▶ If required: Seal the unused connections on the device (1) using the corresponding covers or protective caps.

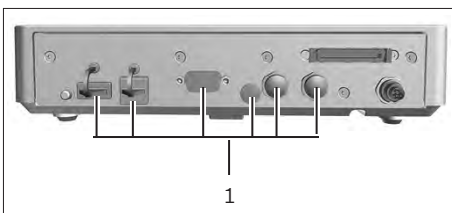

7 System Settings

7.1 Switching the Device On/Off and Activating Standby Mode

If the device is being connected to the power supply for the first time, or after restoring the factory settings: The device is switched on and the Setup Wizard opens. All steps in the Setup Wizard must be completed.

Requirements

The device is connected to the power supply.

Procedure

- ▶ **NOTICE** Pointed or sharp-edged objects may damage the operating display! Only touch the operating display with your fingertips.
- ▶ If the Setup Wizard is shown: Follow the instructions in the Setup Wizard in the operating display.
- ▶ If the login display is opened: Log into the device using a user profile.
- ▶ **NOTICE** Damage to the device when disconnecting the device from the weighing module or electronics module! The device must **not** be disconnected from the weighing module or electronics module during operation.
- ▶ To activate standby mode: Press the [On | Off] button.
- ▶ The [Select standby mode] display is opened.
- ▶ Press the required standby mode.
- ▶ The device displays the time.
- ▶ To switch the device off: There are 2 options.
 - ▶ In the [Select standby mode] display, press the [SWITCH OFF DEVICE] button.
 - ▶ Disconnect the device from the power supply.
- ▶ To switch the device back on again from standby mode or after software-controlled switch-off: Press the on key on the back of the device.
- ▶ To switch the device back on again after disconnection from the power supply: Connect the device to the power supply.




7.2 User Login or Logout

The user selection is only displayed when at least one user is created.

Procedure

- ▶ Press the user selection (1).
- ▶ Select a user, e.g., Administrator.
- ▶ Press the [Login] button.
- ▶ If a password is assigned: The input screen for the password opens.
- ▶ Enter the password and confirm.
- ▶ To log out the active user profile from the device: Press the [Logout] button.
- ▶ If required: Log in a new user.

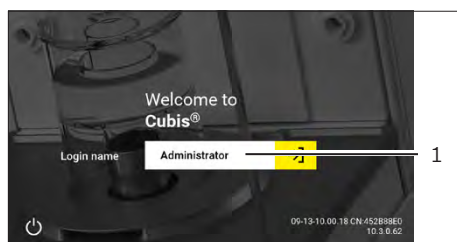

7.3 Performing System Settings

Default settings can be adjusted for the device and the applications in order to align with the ambient conditions and individual operating requirements.

The following settings are necessary to operate the device together with connected components:

- Set up the communication of the connected devices
- Set up additional components

The following settings are recommended to set up the device:

- Set the behavior of the isoCAL function.
- Set the behavior of the motorized draft shield (only for devices with a motorized draft shield).
- If the associated QAPP is activated and the LDAP server is configured: Assign a password.

Procedure

- ▶ Open the main menu.
- ▶ Press the [Settings] button.
- ▶ To adjust settings: Open the preferred submenu.
- ▶ Select the preferred configuration value.
- ▶ Exit the menu.
- ▷ With some settings, the [Booting device] display is displayed in the operating display and the device restarts.

7.4 Using the Help Function

If help texts are available in a menu: The [Help] button is displayed.

Procedure




- ▶ Press the [Help] button.
- ▶ The help texts are displayed.
- ▶ To navigate through the help text: Swipe the text downwards or upwards.

7.5 Activating Applications (QAPPs) and Adding to a Task

7.5.1 Activating Applications

Some applications are activated for the device at the factory. Additional applications may be activated in the QAPP Center. These applications can be tested for 30 days free of charge, and after that require a license.

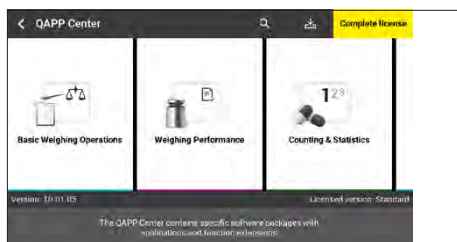

Procedure

- ▶ Open the task.
- ▶ Press the [QAPP Center] button.
- ▷ An overview of the available QAPP packages is displayed.
- ▶ Select the required QAPP package.
- ▷ A list of the applications contained in the QAPP package is displayed.
- ▶ If the selected QAPP package is to be approved with all the applications it contains:
 - ▶ Press the [License] button.
 - ▷ The input field for the license key is displayed.
 - ▶ If an additional cost is associated with the QAPP package: Enter the license key in the input field and press the [OK] button.
 - ▶ If no additional cost is associated with the QAPP package: Press the [OK] button.
- ▶ If an individual application from the displayed QAPP package is to be activated:
 - ▶ Press the preferred application.
 - ▷ A display opens containing details about the selected application.
 - ▶ Press the [License] button.
 - ▷ The input field for the license key is displayed.
 - ▶ If an additional cost is associated with the application: Enter the license key in the input field and press the [OK] button.
 - ▶ If no additional cost is associated with the application: Press the [OK] button.

7.5.2 Adding an Application to a Task

Applications must be added to a task so that they can run.

Procedure

- ▶ Open the task.
- ▶ Press the [New] button.
- ▷ A list of all activated applications is displayed.
- ▶ To select an application: Press the preferred application.
- ▷ The wizard for creating a new task starts.
- ▶ Follow the wizard's instructions in the operating display.

7.6 Switching Off the isoCAL Function

M

If the isoCAL function is switched off for a conformity-assessed device: The device can only be used in restricted temperature ranges for legal-for-trade applications (see Chapter "15.3.2 Ambient Conditions at the Installation Site", page 50). The isoCAL function **cannot** be switched off for all model versions.

Procedure

- ▶ In the "isoCal execution module" submenu, for the "isoCAL function" parameter, select the "Off" configuration value.

7.7 Configuring Motorized Opening and Closing of the Draft Shield (Only for Devices with a Motorized Draft Shield)


Activating the touch sensors on the control unit enables motorized opening and closing of the draft shield. The draft shield has a learning capability. This makes it possible to define how far the draft shield should open when activating the touch sensors. The opening position of the draft shield can be set individually for both touch sensors.

Requirements

The application for use of the motorized draft shield is activated.

Procedure

- ▶ Open the draft shield settings menu.
- ▶ Manually turn the draft shield to the preferred opening angle.
- ▶ Activate the required touch sensor.
- ▷ The draft shield closes.
- ▷ The selected touch sensor is programmed to open the draft shield to the set opening angle.
- ▶ If required: Define a different draft shield opening angle for the other touch sensor.

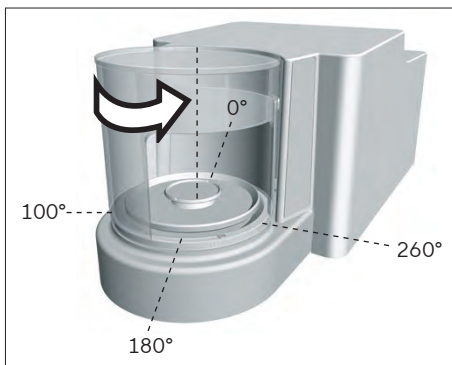

7.8 Adding Weighing and Print Profiles to a Task

To be able to use a weighing or print profile: Add a weighing or print profile to a task. Weighing and print profiles can be configured in the Settings menu.

Procedure

- ▶ Open the task.
- ▶ Create or edit a task. To do this, start the wizard to create or edit a task and follow the wizard's instructions in the operating display.

7.9 Downloading Additional Information

As part of the Cubis® IIICUB firmware package on the Sartorius website, additional information is available for the device, e.g., a description of interface protocols, or a set of installation instructions for a website certificate. The information is available as a PDF file, partially in English. The information is stored on the MySartorius website. A Sartorius ID is required to be able to log into MySartorius and call up information. The Sartorius ID can be generated online.

Procedure

- ▶ Download the “Cubis® CUB Firmware” from the MySartorius website.
- ▶ Find the additional information you require, e.g., the description of interface protocols.

7.10 Observing Warm-up Time

After connecting to the power supply, the warm-up time must be observed. This enables the device to reach its required operating temperature and ensures accurate values during weighing processes.

M

If this relates to a conformity-assessed device: The weight value is marked as **invalid** during the warm-up period.

Procedure

- ▶ Ensure that the warm-up time is observed (see Chapter “15.6 Warm-up Time to Reach Operating Temperature”, page 52).

8 Operation

8.1 Opening and Closing the Motorized Draft Shield (Only for Devices with a Motorized Draft Shield)

8.1.1 Opening or Closing by Touching the Touch Sensors

Requirements

The motorized opening and closing of the draft shield is configured (see Chapter 7.7, page 35).

Procedure

- ▶ Activate the required touch sensor on the right or left side of the control unit. This enables motorized opening or closing of the motorized draft shield according to the stored setting.

8.2 Leveling the Device

Leveling compensates any inclines at the device's installation site. If leveling is necessary: The [Leveling] button is displayed in the weighing display and a message is displayed in the status center.

Procedure

- ▶ If the weighing display is opened: Press the [Leveling] button.
- ▶ When the status center is displayed: Press the [Level] button.
- ▷ The Leveling Wizard opens.
- ▶ Follow the wizard's instructions.

8.3 Calibration, Adjustment or Linearization

Function	Description
Calibration	The device checks how much the displayed value deviated from the specified setpoint.
Adjustment	The device corrects the deviation to the setpoint.
Linearization	The device corrects the deviation from the ideal characteristic curve and the setpoint.

The device needs to be calibrated and adjusted regularly. Various methods can be selected for this:

- Adjusting with the isoCAL function
- Internal calibration and adjustment
- External adjustment
- Internal linearization

M

Only internal adjustment is described below.

External adjustment is **not** possible for conformity-assessed devices in legal metrology.

Procedure

- ▶ Where one of the following conditions occurs, calibrate and adjust the device using the preferred method:
 - Daily, every time the device is switched on
 - After every leveling
 - After changing the ambient conditions (temperature, humidity, or air pressure)
 - After setting the device up at a new installation site

8.3.1 Adjusting With the isoCAL Function

The device can be automatically internally calibrated and adjusted using the isoCAL function.

Requirements

- The isoCAL function is configured in the “Safe weighing” menu, e.g., “On, automatic execution”.
- The conditions for triggering and executing the isoCAL function are met (see Chapter “15.9 Conditions for isoCAL Function”, page 53).

Procedure

- ▶ If the isoCAL function is set to automatic start and the isoCAL function is triggered:
 - ▷ The [isoCAL] button flashes in the operating display.
 - ▶ Wait until the isoCAL function is executed.
 - ▷ In the operating display, a time display counts down from 15 seconds to 0.
 - ▷ If **no** load change or **no** operation takes place on the device before the expiration of the time display: The isoCAL function starts.
- ▶ If the isoCAL function is set to manual start and the isoCAL function is triggered:
 - ▷ The [isoCAL] button flashes in the operating display.
 - ▶ Press the [isoCAL] button.
 - ▷ The isoCAL function starts.
- ▷ If the isoCAL function is complete: The device confirms the completion of the calibration and adjustment process with an acoustic signal, and the calibration report is displayed.
- ▶ To output the calibration report via a connector: Press the [Print Memory] button.
- ▶ To close the calibration report and return to the previous display: Press the [OK] button.

8.3.2 Internally Calibrating and Adjusting the Device**Requirements**

The weighing pan is unloaded.

Procedure

- ▶ Open the main menu.
- ▶ Press the "Internal adjustment" task.
- ▷ The internal calibration and adjustment function is executed.
- ▷ If automatic leveling is configured: The device levels itself automatically.
- ▷ If the calibration and adjustment function is complete: The device confirms the completion of the calibration and adjustment process with an acoustic signal, and the calibration report is displayed.
- ▶ To output the calibration report via a connector: Press the [Print Memory] button.
- ▶ To close the calibration report and return to the previous display: Press the [OK] button.

8.4 Preparing Weighings

The device must be prepared before every weighing.

Procedure

- ▶ Level the device.
- ▶ Zero the device. To do this, press the [Zero] button.
- ▶ If the device **cannot** be zeroed: Remove the sample to be weighed and re-zero the device.
- ▶ Adjust the device.

8.5 Weighing

When weighing chemicals, suitable vessels must be used for the material to be weighed. This makes it possible to prevent damage to the device or its accessories.

Requirements

The device is leveled and adjusted.

Procedure

- ▶ Start a task with weighing function.
- ▶ Zero the device. To do this, press the [Zero] button.
- ▶ If below-balance weighing is being carried out: Suspend the sample on the below-balance weighing hook, e.g., with a wire.
- ▶ If using a vessel for the sample:
 - ▶ Place the vessel for the material to be weighed onto the weighing pan.
 - ▶ Press the [Tare] button. This compensates for the weight of the vessel.
 - ▶ Tare the device. To do this, press the [Tare] button.
 - ▶ Place the sample in the vessel or fill the vessel.
- ▶ If **no** vessel is used for the sample: Place the sample on the weighing pan.
- ▶ Once the weight value is displayed in black and the weighing unit is displayed: Read off the measured value.

8.6 Running Application (Example)

8.6.1 Executing the “Unit change” Function

The “Unit change” function enables switching between the different units and resolutions defined in the weighing profile of the active task. The units and resolutions can be set at the beginning of the weighing process.

Procedure

- ▶ Start the preferred task.
- ▶ Press the [Unit change] button.
- ▷ All units defined in the weighing profile for the active task are displayed in a list.
- ▷ All resolutions for the weight value defined in the weighing profile for the active task are displayed in a list.
- ▶ Press the preferred unit.
- ▶ To set the resolution for the selected unit: Press the preferred resolution.
- ▶ To confirm the selection and return to the weighing display: Press the [OK] button.
- ▷ The current weight value is displayed in the selected unit and resolution.

9 Cleaning and Maintenance

9.1 Preparing the Device for Cleaning

Procedure

- ▶ If an accessory is connected to the device: Disconnect the accessory from the device (see instructions for the accessory).
- ▶ If advanced cleaning needs to be performed: Switch off the device.
- ▶ Open the draft shield.




9.2 Cleaning the Device

Sartorius recommends cleaning the device at regular intervals, e.g., weekly. There must be **no** foreign substances, e.g., particles, fibers or liquids, in the area of the weighing pan.

To clean the device, you can use Sartorius cleaning utensils or a moistened cloth. The device can be cleaned with the help of the Cleaning QAPP. The Cleaning QAPP contains guided processes for normal and advanced device cleaning. The help texts in the Cleaning QAPP contain notes on valid cleaning agents and cleaning intervals (if cleaning intervals have been set).

Procedure

- ▶ **NOTICE** Damage to the device due to an unsuitable tool! The weighing system is hermetically separated from the area of the draft shield plate. This area must **not** be cleaned. This **prevents** dirt from entering.
- ▶ **⚠ WARNING** Risk of injury due to electrical current! Protect the power supply unit and the power supply cable from liquids.
- ▶ Only use suitable cleaning agents and cleaning procedures and observe the product information for the cleaning agent used (cleaning agents see Chapter "15.14 Cleaning Agents and Cleaning Procedures", page 55).
- ▶ If the Cleaning QAPP is to be used for cleaning: Open the Cleaning QAPP for cleaning the device and follow the instructions in the operating display.

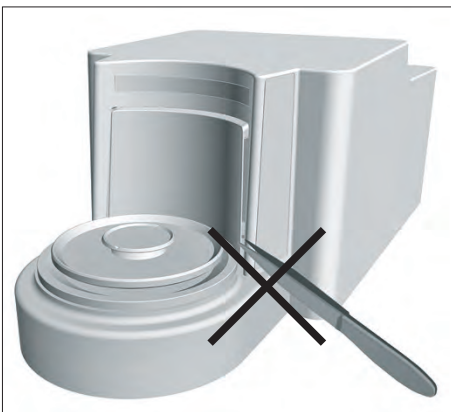

9.3 Recommissioning

Procedure

- ▶ Re-insert all components into the device (see Chapter 5.5, page 26, Chapter 5.6, page 27).
- ▶ Connect the required accessories (see Chapter 6.5, page 31).
- ▶ If the device was disconnected from the power supply: Re-connect the device to the power supply (see Chapter 6.4, page 30).

9.4 Performing a Software Update

A software update can be installed (software package) from a USB mass storage device using one of the device's USB connections. An update can also be carried out from a server via other connectors of the device. The following describes the installation of a USB mass storage device, taken from the Sartorius website.

A software update can extend or change the functionality of the device. Sartorius recommends carrying out software updates:

- Prior to starting the software updates, save the device data to a USB mass storage device.
- If a QAPP Center update is to be carried out in addition: Perform the software update for the device first.

Two files are required for the software update: Firmware file with file ending ".upd", and checksum file with file ending ".upd.md5".

Execution of the software update and troubleshooting during this process are described in the help text "Device maintenance".

Requirements

- The device is connected to the power supply.
- The logged-in user has administrator rights.

Procedure

- ▶ Download the software package from the Sartorius website onto the USB mass storage device. To do this, download the "Cubis® III CUB Firmware" file.
- ▶ If this relates to a zip file: Unzip the software package on the USB mass storage device. The files must be saved in the main directory (at root level). The files must **not** be copied into a folder.
- ▶ Insert the USB mass storage device with the software package into the device's USB-A connection.
- ▶ In the "Settings" / "Device maintenance" menu, press the "Update firmware" menu entry.
- ▶ Press the "USB stick" as the connector and select the preferred software version.
- ▷ The software update takes approx. 3 minutes.
- ▷ Once the software update is complete: The software version number is updated in the login display.

9.5 Performing a QAPP Center Update

The QAPP Center package can be installed from a USB mass storage device using one of the device's USB connections. An update can also be carried out from a server via other connectors of the device. The following describes the installation of a USB mass storage device, taken from the Sartorius website.

Sartorius recommends the following when carrying out the QAPP Center update:

- Prior to starting the QAPP Center update, save the device data to a USB mass storage device.
- If a software update is also to be carried out: Perform the software update for the device first.

Two files are required for the QAPP Center update: QAPP Center with file ending ".appcenter", and checksum file with file ending "qappcenter.upd.md5".

Execution of the QAPP Center update and troubleshooting during this process are described in the help text "Device maintenance".

Requirements

- The device is switched on.
- The QAPP Center package is saved on a USB mass storage device or on a server via a connector.

Procedure

- ▶ Download the QAPP Center package from the Sartorius website onto the USB mass storage device. To do this, download the "Cubis® III CUB Firmware" file.
- ▶ If this relates to a zip file: Unzip the QAPP Center package on the USB mass storage device. The files must be saved in the main directory (at root level). The files must **not** be copied into a folder.
- ▶ Insert the USB mass storage device with the QAPP Center package into one of the device's USB-A connections.
- ▶ In the "Settings" / "Device maintenance" menu, press the "Install QAPP Center" menu entry.
- ▶ Press "USB stick" as the connector.
- ▶ Press the preferred package.
- ▶ Once the QAPP Center update is complete: Confirm successful installation with the [OK] button.
- ▷ Existing tasks remain unchanged after the QAPP Center update. In the existing tasks, the original QAPP versions are used.
- ▶ In order to use the new QAPP version: Create a new task with the new QAPP version. Existing tasks are **not** automatically adapted via the QAPP Center update.

10 Malfunctions

10.1 Status Messages

Status Message	Malfunction	Cause	Solution	Chapter, Page
The LED strip on the weighing module is illuminated in orange.	The device displays an error message.	The device has detected a malfunction.	Rectify the error shown in the operating display.	

10.2 Warning Messages

Warning Message	Malfunction	Cause	Solution	Chapter, Page
Disp.Err.	The value to be output cannot be shown in the operating display.	The data to be displayed is not compatible with the set display format.	Adjust the display settings in the menu, e.g., resolution, unit, decimal places.	
High	The device is overloaded.	The device's maximum weighing capacity has been exceeded.	Reduce the applied weight to below the device's maximum weighing capacity.	15.8, 53
Low	The modulation of the weighing module's weighing converter is too low.	No weighing pan has been placed on the balance. A previously forgotten weight was removed after starting the device.	Insert the weighing pan into the device and switch the device off and on again.	
Com.Err.	The device is not receiving any weight values.	No communication exists between the control unit and the weighing module.	Wait until the control unit restores the communication with the weighing module. If the problem reoccurs: Contact Sartorius Service.	17, 60

10.3 Troubleshooting

Malfunction	Cause	Solution	Chapter, Page
The operating display is black.	The device is disconnected.	Check the connection to the power supply.	6.4, 30
	The power supply unit is not connected.	Connect the power supply cable to the power supply.	6.4, 30
A peripheral device in the USB port is not working.	The connection to the device has been interrupted.	Switch the device off and on again.	7.1, 32
The operating display is red.	The connection to the device has been interrupted.	Switch the device off and on again.	7.1, 32
The displayed weight value changes constantly.	The device installation site is unstable.	Adjust the parameters in the "Ambient conditions" submenu.	
		Change the installation site.	5.2, 24
The weight readout displayed by the device is obviously wrong.	A foreign object is positioned between the weighing pan and the housing.	Remove the foreign object.	
The weight readout displayed by the device is obviously wrong.	The device has not been adjusted.	Adjust the device.	8.3, 37
	The device was not tared before weighing.	Tare the device.	
The motorized draft shield opens or closes in an uncontrolled manner.	An object is permanently in the reception range of the proximity sensor.	Remove the object causing the disturbance or adjust the sensitivity of the proximity sensors in the menu.	

10.4 Operating Issues

Malfunction	Cause	Solution	Chapter, Page
Password forgotten.	A user has forgotten the password.	Contact the administrator to delete or reset the password.	
	The administrator has forgotten the password.	Contact Sartorius Service to delete or reset the password.	17, 60

11 Decommissioning

11.1 Decommissioning the Device

Procedure

- ▶ Disconnect the device from the power supply. To do this, disconnect the power supply cable from the wall outlet.
- ▶ Disconnect all cables and accessories from the device connections.
- ▶ Place all caps and attachment hoods onto the corresponding connections.
- ▶ Clean the device (see Chapter 9.2, page 41).
- ▶ If below-balance weighing has been set up:
 - ▶ Turn the device on its side and place it on a soft surface.
 - ▶ Insert cover plate (1) and screw on with both screws (2).
 - ▶ Place the device back on the device base on a level surface.



12 Transport

12.1 Transporting the Device

Procedure

- ▶ **⚠ CAUTION** Risk of injury from breaking glass! Glass components can break if they fall or are handled incorrectly. Glass fragments can cause cuts.
 - ▶ Only lift the device by its base, **not** by the draft shield.



- ▶ When lifting and transporting, ensure that **no** personnel or objects are in the way.
- ▶ Get help from others when lifting and transporting the device.
- ▶ Use suitable conveyance devices, e.g., trolleys, for long transport routes.

13 Storage and Shipping

13.1 Storing

Procedure

- ▶ Decommission the device.
- ▶ Store the device according to the ambient conditions (see Chapter “15.2 Ambient Conditions During Storage and Transport”, page 49).

13.2 Returning the Device and Parts

Defective devices or parts can be returned to Sartorius. Returned devices must be clean and packed in their original packaging.

Transport damage as well as measures for subsequent cleaning and disinfection of the device or parts by Sartorius shall be charged to sender.

Devices contaminated with hazardous materials, e.g., harmful biological or chemical substances, will **not** be accepted for repair or disposal.

Procedure

- ▶ Decommission the device.
- ▶ Contact Sartorius Service for instructions on how to return devices or parts.
- ▶ Pack the device and its parts in their original packaging for return.

14 Disposal

14.1 Disposing of the Device and Parts

The device and the device accessories must be disposed of properly by disposal facilities.

Two lithium cell batteries, type CR2032, are installed inside the device. Batteries must be disposed of properly by disposal facilities.

Many of the packaging materials are recyclable to promote sustainability and to contribute to reducing global waste volumes.

Procedure

- ▶ Dispose of the device in accordance with local government regulations. Inform the disposal facility that there are two lithium cell batteries, type CR2032, installed inside the device.
- ▶ Dispose of the packaging in accordance with local government regulations.

15 Technical Data

15.1 Dimensions and Weight

	Unit	With round glass ring draft shield	With stainless steel draft shield for the filter weighing pan
		Value	Value
Dimensions			
Weighing module (L × W × H)	mm	340 × 139 × 129	340 × 139 × 129
Electronics module with control unit (L × W × H)	mm	340 × 240 × 61	340 × 240 × 61
Weighing pan size			
Weighing pan for model CUB2.7S	mm	Ø 20	Ø 20
Weighing pan for models CUB10.6S, CUB6.6S, CUB3.6P	mm	Ø 30	Ø 30
Filter weighing pan	mm	–	Ø 50
Weight, approx.	kg	7.2	7.9

15.2 Ambient Conditions During Storage and Transport

	Unit	Value
Temperature		
During storage and transport	°C	-20 – +60
Storage, dry		

15.3 Installation Conditions

15.3.1 Installation Site

	Unit	Value
Height above sea level, maximum	m	3000
No potentially explosive atmospheres		
Laboratory room, with pollution level according to DIN EN 61010-1		2
Suitable for protection class		
Protection class of the device, according to EN 60529-1		IP30
Protection class of the power supply unit, according to DIN EN 60529		IP30
Access to operation-relevant parts is guaranteed		

	Unit	Value
Space requirements		
Suitable for the dimensions of the device and the associated components		
Setup surface		
Stable, even		
Suitable for the weight of the device and the associated components		
Additional properties		
No heat from heating systems or direct sunlight		
No direct drafts from open windows, AC systems, or doors		
No vibrations		
No "heavy traffic" areas (personnel)		
No electromagnetic fields or electromagnetic radiation, e.g., from radio equipment		

15.3.2 Ambient Conditions at the Installation Site

	Unit	Value
Temperature		
In operation	°C	+5 - +40
In operation, with isoCAL function, scope of application as per Directive 2014/31/EU	°C	+10 - +30
In operation, without isoCAL function, scope of application as per Directive 2014/31/EU	°C	+17 - +27
In operation, with conformity-assessed devices, as per specifications on the device ID label		
Relative humidity, in operation		
At temperatures up to 31°C, maximum	%	80
Then linear decrease from 80 % at 31 °C to 50 % at 40 °C		

15.4 Acclimatization Before Power Supply

	Unit	Value
Time period between unpacking and connecting to a power supply	h	2

15.5 Electrical Data

15.5.1 Power Supply

	Unit	Value
Power supply only permitted via the power supply cable and power supply unit supplied by Sartorius		
Sartorius power supply unit, type YEPS03-15V0		
Primary (power supply unit)		
AC voltage	V	100 - 240 ± 10%
Frequency	Hz	50 - 60 ± 5%
Current consumption, maximum	A	1.0
Power consumption, maximum	W	
Secondary (device)		
DC voltage	V	14.25 - 15.75
Current consumption, maximum	A	2
Power consumption, maximum	W	30
Device fuses		
Quantity		1
Type: Electronic		
Power supply unit fuses		
Quantity		1
Type: Electronic		
Protection class, according to IEC 62368-1		
Device		1
Power supply unit		1
Overvoltage category according to IEC 61010-1		
Device		II
Power supply unit		II
Potential equalization, via potential equalization cable, not supplied with the device		
Other data: See label on the power supply unit		

15.5.2 Safety of Electrical Equipment

According to EN 61010-1 / IEC 61010-1 Safety requirements for electrical equipment for measurement, control, and laboratory use – Part 1: General requirements

15.5.3 Electromagnetic Compatibility

Interference immunity
Suitable for use in industrial areas
Transient emissions
Class B
Suitable for use in residential areas and areas that are directly connected to a low voltage network that also supplies residential buildings.

15.6 Warm-up Time to Reach Operating Temperature


	Unit	Value
Time period between switching the device on and carrying out weighing tasks	h	1

15.7 Interfaces

15.7.1 Specifications for the COM-RS232 Interface

Type of interface	Serial interface
Interface operation	Full duplex
Level	RS232
Connection	D-sub connector, 9-pin

Pin assignment

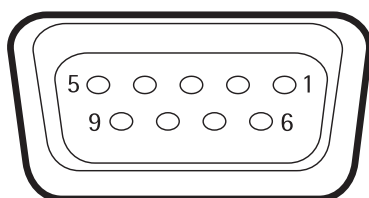

- Pin 1: **Not** assigned
- Pin 2: Data output (TxD)
- Pin 3: Data input (RxD)
- Pin 4: **Not** assigned
- Pin 5: Internal ground
- Pin 6: **Not** assigned
- Pin 7: Clear to Send (CTS)
- Pin 8: Request to Send (RTS)
- Pin 9: **Not** assigned

15.7.2 Specifications for the USB-A Interface

Communication	USB host (Master)
Connectable devices	Sartorius printers, USB memory sticks, USB barcode scanners, USB keyboards, NFC readers

15.7.3 Specifications for the USB-B Interface

Communication	USB device (Slave)
Type of interface	Virtual serial interface (virtual COM-port, VCP) and "PC-Direct" communication

15.7.4 Specifications for the USB-C Interface

Communication	Downstream-facing port (DFP), USB host (Master)
Communication	RS232 connection with accessory YCC-USB-C-D09M

15.8 Recommended Calibration Weight

		CUB2.7S	CUB10.6S	CUB6.6S	CUB3.6P
	Unit	Value	Value	Value	Value
External test weight	g	2	10	5	3
Recommended accuracy class		E2	E2	E2	E2

15.9 Conditions for isoCAL Function

	Unit	Value
Possible conditions for triggering the isoCAL function		
In the event of a temperature change	K	1.5
After a time interval	h	12
After successful leveling		
Only conformity-assessed devices: After interrupting the power supply		
Required conditions for executing the isoCAL function		
Operating display is in weighing mode (not in menu)		
Alphanumeric inputs are not active		
Min. time period without entry on the device	min	2
Min. time period with unchanged loading of the weighing pan	min	2
Max. loading of weighing pan, high-capacity	%	2

15.10 Data Storage Device

	Value
Maximum number of entities	500000

15.11 Integrated Clock

	Unit	Value
Maximum deviation per month (RTC)	s	30

15.12 Backup Battery

	Unit	Value
Lithium battery, type CR2032		
Service life at room temperature, minimum	Years	10

15.13 Materials

Housing:

Die-cast aluminum

Plastic PBT

Float glass Optiwhite

Control unit

Plastic PA12

Operating display: Plastic PA12 | float glass

15.14 Cleaning Agents and Cleaning Procedures

15.14.1 Approved Cleaning Agents

Device components	Cleaning agents and concentration					
	Ethanol, 70%	Isopropanol, 70 %	Citric acid, 10 %	Diluted hydrogen peroxide, 3.5%	Sodium hydroxide, 32%	Ecolab Klercide™ Sporicidal Chlorine Thiosulfate
Draft shield	x	x	xx	xx	-	xx
Components in the weighing compartment						
Filter weighing pan	x	x	x	x	xx	x
Weighing pan	x	x	x	x	xx	x
Shield plate	x	x	x	x	xx	x
Rear wall of weighing compartment	xx	x	x	x	x	x
Control unit	x	x	x	x	x	x
Back of the device						
Plastic surfaces	x	xx	x	x	x	x
Electronics module	x	xx	x	x	x	x
x Suitable						
xx Suitable, may affect visual appearance of the device, does not impact mechanical stability						
- Not suitable						

15.14.2 Approved Cleaning Procedures

Remove dust and powdery sample residue with a brush or hand-held vacuum cleaner.

Use blotting paper to remove any liquid sample material.

Wipe the components of the weighing pan with a brush or a lightly moistened cleaning cloth.

Spraying the device surfaces with cleaning agents, contact time Min 5-10

Wiping the device surfaces with a lightly moistened cleaning cloth; for more severe contaminations, use a mild soap solution.

15.15 Metrological Data

15.15.1 Models CUB2.7S | CUB10.6S | CUB6.6S | CUB3.6P

		CUB2.7S	CUB10.6S	CUB6.6S	CUB3.6P
	Unit	Value	Value	Value	Value
Scale interval (d)	µg	0.1	1	1	1 2 5
Maximum capacity (Max)	g	2.1	10.1	6.1	1.1 2.1 3.1
Repeatability at 5 % load					
Standard deviation of the load values, tolerance	µg	0.20	1	1	3
Standard deviation of the load values, typical value	µg	0.15	0.4	0.4	0.4
Repeatability at approx. maximum capacity					
Standard deviation of the load values, tolerance	µg	0.25	1	1	5
Standard deviation of the load values, typical value	µg	0.18	0.6	0.6	0.6
Linearity deviation					
Tolerance	µg	0.9	4	4	4
Typical value	µg	0.7	3	3	3
Deviation when load is off-center, positions according to OIML R76					
Test weight	g	1	5	2	1
Tolerance	µg	0.7	4	4	5
Typical value	µg	0.5	3	3	3
Sensitivity drift between +10 °C and +30 °C	ppm/K	1	1	1	1
Tare-maximum capacity: Less than 100 % of maximum capacity					
Accuracy class, according to Directive 2014/31/EU		I	I	I	I
Verification scale interval (e), according to Directive 2014/31/EU	mg	1	1	1	1
Minimum load (Min), according to Directive 2014/31/EU	mg	0.01	0.1	0.1	0.1
Minimum initial weight according to USP (United States Pharmacopeia), Chap. 41					
Optimum minimum initial weight	µg	82	820	820	820
Typical minimum initial weight	µg	300	820	820	820
Typical stabilization time	s	7	3	3	3
Typical measurement time	s	10	6	6	6

16 Accessories

16.1 Accessories

This table contains a selection of the accessories that can be ordered. For information on other products, contact Sartorius.

16.1.1 Printers and Communication

Item	Quantity	Order Number
Thermal transfer direct thermal printer for GLP GMP printouts on continuous paper and tags	1	YDP30
Thermal transfer direct thermal network printer with Ethernet connection for GLP GMP printouts on continuous paper and tags	1	YDP30-NET
Printer paper and thermal transfer ribbon, 90 m	1	69Y03285
Self-adhesive paper and thermal transfer ribbon, 90 m	1	69Y03286
Thermal paper, 24 m	5	69Y03287
Self-adhesive thermal paper, 24 m	5	69Y03288
Tag roll, 58 x 30 mm	1000	69Y03092
Tag roll, 58 x 76 mm	500	69Y03093
Tag roll, 58 x 100 mm	350	69Y03094
Wireless nano USB adapter for company network or independent Wi-Fi network, e.g., operation with a Sartorius network printer YDP30-NET (only for Europe)	1	YWLAN01MS
Wireless nano router, e.g., for Sartorius network printer YDP30-NET for operation in an independent Wi-Fi network (only for Europe)	1	YWLAN02MS
Display cable, 3 m, for separate installation of display and weighing units, installation by Sartorius Service or at the factory	1	YCC01-CUB-2
Installation of display cable, 3 m, for separate installation of display and weighing units	1	VF4016
Cable, 3 m, between the weighing module and electronics module for semi-microbalance	1	YCC01-MSM3
USB cable, 3 m, for connecting USB B port to USB A port	1	69MS0099
RS232C connection cable, 9-pin, 3 m, for connection to a PC with 9-pin COM interface	1	VF4761
USB QR barcode scanner	1	YBR05
USB-RFID reader	1	YRFID01
Foot switch for the "open closed" draft shield functions (only in combination with motorized draft shield), taring and printing	1	YFS02
Motion sensor for triggering a maximum of 4 functions via gesture control, selection via menu	1	YHS02USB
Dust cover	1	YDCC2CUB

16.1.2 Hardware and Software for Pipette Calibration

Item	Quantity	Order Number
Pipette calibration kit (hardware) for microbalance weighing modules, consists of moisture trap and all required adapters	1	YCP08MC

16.1.3 Titanium Holder

Item	Quantity	Order Number
Titanium holder		
For coronary stents up to 38 mm	1	YSH10
For safe-lock tubes, 1.5 – 2 ml	1	YSH13

16.1.4 Filter Balance and Antistatic Accessories

Item	Quantity	Order Number
Titanium filter weighing pan, 52 mm diameter, for filter microbalances, ultramicrobalances or microbalances	1	YSH34
Titanium filter weighing pan, 75 mm diameter, for filter microbalances, ultramicrobalances or microbalances	1	YSH35
Titanium filter weighing pan, 90 mm diameter, for filter microbalances, ultramicrobalances or microbalances	1	YSH36
Ionizer, including universal adapter	1	YIB03-C

16.1.5 Special Applications

Item	Quantity	Order Number
Cleaning kit	1	YCK01MC

16.1.6 Weighing Tables

Item	Quantity	Order Number
Weighing table		
Made from synthetic stone, with vibration dampening	1	YWT03
Made from wood and synthetic stone	1	YWT09
Wall console	1	YWT04

16.1.7 Weighing Accessories

Item	Quantity	Order Number
Aluminum weighing scoop, 4.5 mg for ultramicrobalances and microbalances	250	6565-250
Aluminum weighing scoop, 52 mg for ultramicrobalances and microbalances	50	6566-50

17 Sartorius Service

Sartorius Service is available for queries regarding the device. Visit the Sartorius website for information about the service addresses, services provided, or to contact a local representative.

For inquiries about the system or when contacting Sartorius Service in the event of a malfunction, keep the device information close at hand, e.g., serial number, hardware, firmware, configuration, to pass it on to Sartorius Service. To do this, refer to the information on the manufacturer's ID label and in the "General device Information" menu.

18 Trademark Information

Ecolab Klercide™ is a registered trademark of Ecolab Europe GmbH.

19 Conformity

19.1 EU Declaration of Conformity

The attached Declaration of Conformity hereby confirms compliance of the device with the directives cited.

M

The Declaration of Conformity supplied with the balance is for conformity-assessed (verified) balances for use in the EEA. Keep it in a safe place.