




## Personal floor scale KERN MPE • MPE-E

High quality personal floor scale with BMI function – with EC type approval and approval for professional medical use in medical diagnostics

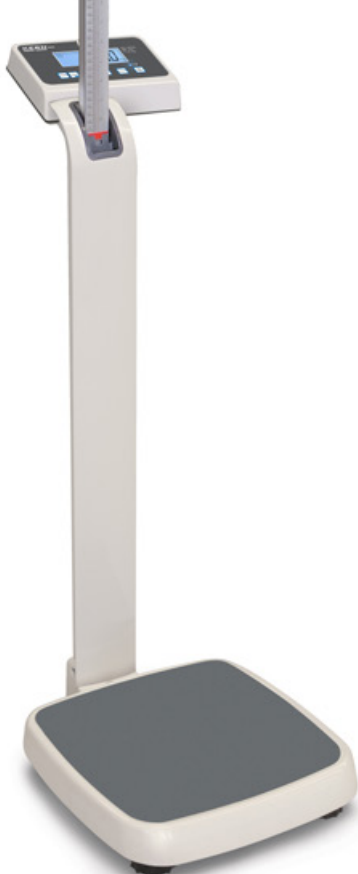

**KERN MPE-H**  
with stand and  
height rod



**KERN MPE-P**  
with stand

## Personal floor scale KERN MPE · MPE-E



### Features

- Verification class III (verification is optional)
- Approved as a medical device according to 93/42/EEC
- **1** Practical second display on the back of the scale, which makes it easy for both the patient and the doctor or nursing staff to read the weight
- **2** The MPE-H model has an elegant height rod integrated into the column, which can be pulled out for use and has a sturdy plastic head rest and a measuring range from 88-205 cm
- The ergonomically optimised shape and robust design mean that this range is ideal for daily use in a professional environment
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely
- Large platform with non-slip and wear-resistant surface made of plastic, to provide a secure footing for patients
- Hold function: When patients do not stand completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Mother and child function: the parent's weight will be taken and then the scale will be reset to zero at the push of a button. Afterwards then the parent and baby can be weighed together. The weight of the child will now be shown on the display. This is particularly good for infants, as this way they can be cradled, safe and sound, in the arms of the parent while being weighed

- BMI function to determine underweight/normal weight/surplus weight
- High resolution readability: Readability [d] can be increased by one decimal place for 5 sec. by the touch of a key
- Data interface RS-232 as standard
- **3** Easy to transport thanks to the two integrated transport wheels
- Battery- or mains-powered, rechargeable battery operation optional
- Protective working cover included with delivery
- KERN MPE-E: Identical in construction to Kern MPE, however without mains adapter, only battery operation. With a particularly attractive price

### Technical data

- Large backlit LCD display, digit height 25 mm
- Weighing plate dimensions W×D×H 365×360×80 mm
- Overall dimensions W×D×H MPE-PM: 365×570×1030 mm MPE-HM: 365×570×2134 mm, with extracted height rod
- Net weight approx. 12 kg

### KERN MPE

- Battery operation possible, 6×1.5 V AA not included, operating time up to 20 h
- Mains adapter external, standard

### KERN MPE-E

- Batteries included, 6×1.5 V AA, operating time up to 20 h


### Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MPE-A01S05
- Cleaning cloths, alcohol-free cloths for disinfectant wiping, quick acting, based on modern quaternary ammonium compounds and effective against papova viruses. Particularly gentle on materials, and very well suited for disinfecting products which are sensitive to alcohol. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 100 pcs., size 20×22 cm per cloth, KERN MYC-01
- Internal rechargeable battery pack, operating time up to 40 h without backlight, charging time approx. 12 h, KERN MBC-A08
- External mains adapter, 100 V - 240 V, Standard EU, CH, UK, KERN YKA-43
- External mains adapter, 100 V - 240 V, Standard EU, KERN YKA-44
- Interface cable RS-232 to connect an external device, KERN CFS-A01
- Matrix needle printer, KERN YKN-01
- Universal label printer, KERN YKE-01
- Thermal printer, KERN YKB-01N
- For further details, plenty of further accessories and suitable printers, see Internet

#### STANDARD



#### OPTION

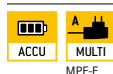

#### FACTORY



Model	Weighing capacity	Readability	Verification value	Mandatory by law Verification
	[Max] kg	[d] kg	[e] kg	<b>KERN</b>
KERN MPE 200K-1PEM	250	0,1	0,1	965-129
MPE 250K100PM	250	0,1	0,1	965-129
MPE 200K-1HEM	250	0,1	0,1	965-129
MPE 250K100HM	250	0,1	0,1	965-129

\* Within the EU, official verification is mandatory by law for scales that are intended for use as a medical device. Please add this to your order. We require the location of use and the post code for the verification.



### Adjusting program CAL:

For quick setting up of the balance's accuracy. External adjusting weight required



### Memory:

Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



### Data interface RS-232:

To connect the balance to a printer, PC or network



### RS-485 data interface:

To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible



### Bluetooth\* data interface:

To transfer data from the balance to a printer, PC or other peripherals



### WIFI data interface:

To transfer data from the balance to a printer, PC or other peripherals



### Control outputs

(optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



### Statistics:

using the saved values, the device calculates statistical data, such as average value, standard deviation etc.



### PC Software:

to transfer the measurements from the device to a PC



### GLP/ISO log:

With date and time. Only with KERN printers



### KERN Communication Protocol (KCP):

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



### Piece counting:

Reference quantities selectable. Display can be switched from piece to weight



### Totalising level A:

The weights of similar items can be added together and the total can be printed out



### Weighing units:

Can be switched to e.g. nonmetric units. Please refer to website for more details



### Weighing with tolerance range

(Check weighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model



### ZERO:

Resets the display to "0"



### Hold function:

When patients do not stand, sit or lie completely still, a stable weight is calculated using an average weight



### Hold function:

When the weighing conditions are unstable, a stable weight is calculated as an average value



### Protection against dust and water splashes IPxx:

The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013



### Suspended weighing:

Load support with hook on the underside of the balance



### Battery operation:

Ready for battery operation. The battery type is specified for each device



### Rechargeable battery pack:

Rechargeable set



### Battery operation rechargeable

Prepared for a rechargeable battery operation



### Universal plug-in power supply:

with universal input and optional input socket adapters for  
A) EU, CH  
B) EU, CH, GB, USA



### Plug-in power supply:

230V/50Hz in standard version for EU. On request GB, AUS or USA version available



### Integrated power supply unit:

Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request



### Weighing principle: Strain gauges

Electrical resistor on an elastic deforming body



### Peak hold function:

capturing a peak value within a measuring process



### Push and Pull:

the measuring device can capture tension and compression forces



### Integrated scale:

In the eyepiece



### 360° rotatable microscope head



### Monocular Microscope:

For the inspection with one eye



### Binocular Microscope:

For the inspection with both eyes



### Trinocular Microscope:

For the inspection with both eyes and the additional option for the connection of a camera



### Abbe Condenser:

With high numerical aperture for the concentration and the focusing of light



### Halogen illumination:

For pictures bright and rich in contrast



### LED illumination:

Cold, energy-saving and especially long-life illumination



### Fluorescence illumination for compound microscopes:

With 100 W mercury lamp and filter



### Fluorescence illumination for compound microscopes:

With 3 W LED illumination and filter



### Phase contrast unit:

For a higher contrast



### Darkfield condenser/unit:

For a higher contrast due to indirect illumination



### Polarising unit:

To polarise the light



### Infinity system:

Infinity corrected optical system



### Automatic temperature compensation:

For measurements between 10 °C and 30 °C



### Verification possible:

The time required for verification is specified in the pictogram



### Package shipment:

The time required for internal shipping preparations is shown in days in the pictogram



### Pallet shipment:

The time required for internal shipping preparations is shown in days in the pictogram

\*The Bluetooth\* word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.