

Polarimeter KERN OAB-L



The ideal helper for getting started with the analysis of your optically active solutions in the laboratory


### Features

- The KERN OAB 10LN is a manual polarimeter which is characterised by its ergonomic design and easy handling
- The powerful 589 nm sodium vapour lamp is the optimum light source to produce a linear, polarised beam of light
- The 1° scale division including Nonius (0.05°) enables precise definition of the angle of rotation of the substance to be observed
- To hold liquid samples, two glass cuvettes (100 mm/200 mm) are included with the delivery
- Included with delivery:  
Sodium vapour lamp, 100 mm glass cuvette, 200 mm Glasküvette, Replacement lenses and sealing rings for cuvettes

### Technical data

- Light source: Sodium vapour lamp (589 nm)
- Stabilisation time: 10 mins after switching on
- Overall dimensions W×D×H  
500×135×330 mm
- Net weight approx. 5 kg

STANDARD

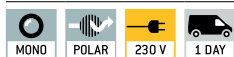

### Polarimeter KERN OAB-L

#### Scope of application: Laboratory/Education

The reliable polarimeters in the OAB-L range have been designed for simple laboratory applications as well as practical training. You can evaluate liquid, optically-active samples with chiral characteristics with this device. Typical applications are determining kinetics in cane sugar inversion, determining mutarotation of glucose and investigation of starch hydrolysis. The optical rotation is measured in degrees.

The main scope of applications is:

- Pharmacy
- Sugar industry: for example cane sugar
- Beverage industry
- Food industry
- Chemical industry
- Laboratories
- Training



Cuvette in measuring chamber




























Model	Scales	Measuring range	Division	Vernier	Wave length
<b>KERN</b>					
<b>OAB 10LN</b>	Optical rotation	± 180°	1°	0,05°	589 nm

#### Accessory parts: OAB

Model	Description
<b>KERN</b>	
<b>OAB-A2501</b>	Glass cuvette, Length: 100 mm (Spare part)
<b>OAB-A2502</b>	Glass cuvette, Length: 200 mm (Spare part)
<b>OAB-A2581</b>	Sodium vapour lamp, Wave length: 589 mm (Spare part)



Cuvette 10 and 20 cm

- 
**360° rotatable microscope head**
- 
**Monocular Microscope**  
 For the inspection with one eye
- 
**Binocular Microscope**  
 For the inspection with both eyes
- 
**Trinocular Microscope**  
 For the inspection with both eyes and the additional option for the connection of a camera
- 
**Abbe Condenser**  
 With high numerical aperture for the concentration and the focusing of light
- 
**Halogen illumination**  
 For pictures bright and rich in contrast
- 
**LED illumination**  
 Cold, energy-saving and especially long-life illumination
- 
**Incident illumination**  
 For non-transparent objects
- 
**Transmitting illumination**  
 For transparent objects
- 
**Fluorescence illumination**  
 For stereomicroscopes
- 
**Fluorescence illumination for compound microscopes**  
 With 100 W mercury lamp and filter
- 
**Fluorescence illumination for compound microscopes**  
 With 3 W LED illumination and filter
- 
**Phase contrast unit**  
 For a higher contrast
- 
**Darkfield condenser/unit**  
 For a higher contrast due to indirect illumination
- 
**Polarising unit**  
 To polarise the light
- 
**Infinity system**  
 Infinity corrected optical system
- 
**Zoom magnification**  
 For stereomicroscopes
- 
**Auto-focus**  
 For automatic control of the focus level
- 
**Parallel optical system**  
 For stereomicroscopes, enables fatigue-proof working
- 
**Integrated scale**  
 In the eyepiece
- 
**SD card**  
 For data storage
- 
**USB 2.0 digital camera**  
 For direct transmitting of the picture to a PC
- 
**USB 3.0 digital camera**  
 For direct transmitting of the picture to a PC
- 
**WIFI data interface:**  
 For transmitting of the picture to a mobile display device
- 
**HDMI digital camera**  
 For direct transmitting of the picture to a display device
- 
**PC software**  
 To transfer the measurements from the device to a PC.
- 
**Automatic temperature compensation**  
 For measurements between 10 °C and 30 °C
- 
**Protection against dust and water splashes IPxx:**  
 The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013
- 
**Battery operation**  
 Ready for battery operation. The battery type is specified for each device.
- 
**Battery operation rechargeable**  
 Prepared for a rechargeable battery operation
- 
**Plug-in power supply**  
 230V/50Hz in standard version for EU. On request GB, AUS or USA version.
- 
**Integrated power supply unit**  
 Integrated in microscope. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.
- 
**Package shipment**  
 The time required to manufacture the product internally is shown in days in the pictogram.

## ABBREVIATIONS

- C-Mount** Adapter for the connection of a camera to a trinocular microscope
- FPS** Frames per second
- H(S)WF** High (Super) Wide Field (Eyepiece with high eye point for wearers of glasses)
- LWD** Long Working Distance
- N.A.** Numerical Aperture
- SLR camera** Single-Lens Reflex camera
- SWF** Super Wide Field (Field number at least  $\varnothing$  23 mm for 10 $\times$  eyepiece)
- W.D.** Working Distance
- WF** Wide Field (Field number up to  $\varnothing$  22 mm for 10 $\times$  eyepiece)