


Tablet Camera KERN ODC-2



ODC 241

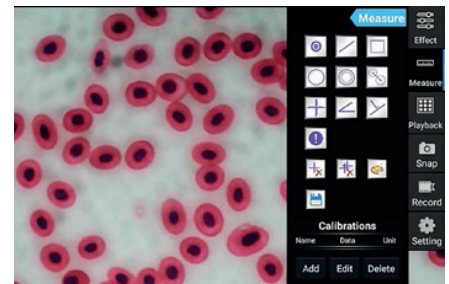

Integrated software with measuring function

Digital microscopy brought up to date – tablet with integrated camera for optimal observation and digital documentation of the sample

Features

- A 2-in-1 solution in digital microscopy as a universal system for trinocular microscopes with C-mount adapter. The ODC 241 microscope-tablet-camera consists of a large Android tablet in combination with a 5-MP camera
- The KERN ODC 241 tablet-camera has been specially developed for simple and direct observation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- As well a live transfer of the image to the Android table, the integrated 5-MP camera also means that images and videos can be created for the documentation.
- Simple measuring functions such as, for example, functions for measuring distance, surfaces and angles as well as a manual counting function are also available
- Automatic white balance and automatic contrast adjustment can be performed quickly and easily, which enables efficient working procedures
- A range of additional functions are provided through the integrated interfaces, such as, for example
 - Data storage on a USB stick or SD card
 - Connection to a USB mouse
 - Transfer of the live image to an external screen using HDMI
 - Transfer of stored data to external receivers using WiFi
- The delivery includes the tablet camera with pre-installed software as well as the mains adapter

Technical data































- 9.7" LCD-Touchscreen
- Screen resolution: 2048×1536 pixels
- CPU: Quad Core Cortex-A17; 1,8 GHz
- Overall dimensions W×D×H 238×51×206 mm
- Net weight approx. 0,65 kg

! Cannot be combined with the microscopes of the OZM-5 series.

STANDARD



Model	Resolution Camera	Interface	FPS	Sensor	Sensor size	Supported operating system
KERN ODC 241	5 MP	WiFi, USB 2.0, HDMI, SD	15 - 30	CMOS	1/2,5"	Android 5.1

- 
360° rotatable microscope head
- 
Monocular Microscope
 For the inspection with one eye
- 
Binocular Microscope
 For the inspection with both eyes
- 
Trinocular Microscope
 For the inspection with both eyes and the additional option for the connection of a camera
- 
Abbe Condenser
 With high numerical aperture for the concentration and the focusing of light
- 
Halogen illumination
 For pictures bright and rich in contrast
- 
LED illumination
 Cold, energy-saving and especially long-life illumination
- 
Incident illumination
 For non-transparent objects
- 
Transmitting illumination
 For transparent objects
- 
Fluorescence illumination
 For stereomicroscopes
- 
Fluorescence illumination for compound microscopes
 With 100 W mercury lamp and filter
- 
Fluorescence illumination for compound microscopes
 With 3 W LED illumination and filter
- 
Phase contrast unit
 For a higher contrast
- 
Darkfield condenser/unit
 For a higher contrast due to indirect illumination
- 
Polarising unit
 To polarise the light
- 
Infinity system
 Infinity corrected optical system
- 
Zoom magnification
 For stereomicroscopes
- 
Auto-focus
 For automatic control of the focus level
- 
Parallel optical system
 For stereomicroscopes, enables fatigue-proof working
- 
Integrated scale
 In the eyepiece
- 
SD card
 For data storage
- 
USB 2.0 digital camera
 For direct transmitting of the picture to a PC
- 
USB 3.0 digital camera
 For direct transmitting of the picture to a PC
- 
WiFi data interface:
 For transmitting of the picture to a mobile display device
- 
HDMI digital camera
 For direct transmitting of the picture to a display device
- 
PC software
 To transfer the measurements from the device to a PC.
- 
Automatic temperature compensation
 For measurements between 10 °C and 30 °C
- 
Protection against dust and water splashes IPxx:
 The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013
- 
Battery operation
 Ready for battery operation. The battery type is specified for each device.
- 
Battery operation rechargeable
 Prepared for a rechargeable battery operation
- 
Plug-in power supply
 230V/50Hz in standard version for EU. On request GB, AUS or USA version.
- 
Integrated power supply unit
 Integrated in microscope. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.
- 
Package shipment
 The time required to manufacture the product internally is shown in days in the pictogram.

ABBREVIATIONS

- C-Mount** Adapter for the connection of a camera to a trinocular microscope
- FPS** Frames per second
- H(S)WF** High (Super) Wide Field (Eyepiece with high eye point for wearers of glasses)
- LWD** Long Working Distance
- N.A.** Numerical Aperture
- SLR camera** Single-Lens Reflex camera
- SWF** Super Wide Field (Field number at least \varnothing 23 mm for 10 \times eyepiece)
- W.D.** Working Distance
- WF** Wide Field (Field number up to \varnothing 22 mm for 10 \times eyepiece)