

## Video microscope KERN OIV-2



OIV 254 Snapshot button

The comprehensive digital solution for increased working comfort when carrying out continuous monitoring work in industry.

### Features Technical data

- The Kern OIV-2 is a video microscope which has been constructed to optimise digital stereo microscopy. Our well-conceived, comprehensive solution with axial optical unit enables immediate, simple display of your samples on the screen.
- The LED incident illumination unit (ring) included as standard guarantees the very best illumination of your sample.
- Combined with the large working surface, recording objects on the screen is ideally suited for monitoring, analysis and documentation in industrial environments.
- The excellent optical unit enables continuous sharp image tracking across the entire zoom range from 0.7×–5×.


- The powerful 2.0 megapixel camera of the microscope without eyepieces offers, thanks to the HDMI output, smooth live monitoring of your samples from the HD monitor. In addition, the software which is easy to use, the USB stick as well as the USB mouse which are integral components of the delivery, mean you can process and store your results digitally.
- With the OIV 254 model, there is the option of image capture at the push of a button, without having to detour via the software. Whereas the OIV 255 guarantees software-controlled taking of images and videos with additional, documentation functions
- A protective dust cover, as well as multi-lingual user instructions are included in the scope of the delivery

- Optical system: Axial
- Brightness adjustable
- Screen: 12", 1920×1080 HD, -5°–15° inclination
- Magnification ratio: 7,1:1
- Stand: arm curved
- Illumination: 2 W LED ring (incident)
- Data storage: External using USB (Max 128 GB)
- Working distance: 105 mm
- Maximum sample height: 100 mm
- Overall dimensions W×D×H 320×260×483 mm
- Net weight approx. 6 kg







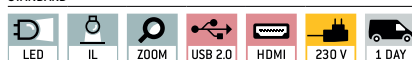























### Accessories

- Auxiliary objective 0,5×, KERN OZB-A2101

STANDARD



Model	Standard configuration					
	Resolution camera	Interface	Sensor	Field of view mm	Objective Zoom	Software functions
<b>KERN OIV 254</b>	2 MP	HDMI (60 FPS)	CMOS 1/2"	∅ 29,82–4,18	0,7×–5×	Image capture
<b>OIV 255</b>	2 MP	HDMI (60 FPS)	CMOS 1/2"	∅ 29,82–4,18	0,7×–5×	Images and videos, documentation

- 
**360° rotatable microscope head**
- 
**Monocular Microscope**  
 For the inspection with one eye
- 
**Binocular Microscope**  
 For the inspection with both eyes
- 
**Trinocular Microscope**  
 For the inspection with both eyes and the additional option for the connection of a camera
- 
**Abbe Condenser**  
 With high numerical aperture for the concentration and the focusing of light
- 
**Halogen illumination**  
 For pictures bright and rich in contrast
- 
**LED illumination**  
 Cold, energy-saving and especially long-life illumination
- 
**Incident illumination**  
 For non-transparent objects
- 
**Transmitting illumination**  
 For transparent objects
- 
**Fluorescence illumination**  
 For stereomicroscopes
- 
**Fluorescence illumination for compound microscopes**  
 With 100 W mercury lamp and filter
- 
**Fluorescence illumination for compound microscopes**  
 With 3 W LED illumination and filter
- 
**Phase contrast unit**  
 For a higher contrast
- 
**Darkfield condenser/unit**  
 For a higher contrast due to indirect illumination
- 
**Polarising unit**  
 To polarise the light
- 
**Infinity system**  
 Infinity corrected optical system
- 
**Zoom magnification**  
 For stereomicroscopes
- 
**Auto-focus**  
 For automatic control of the focus level
- 
**Parallel optical system**  
 For stereomicroscopes, enables fatigue-proof working
- 
**Integrated scale**  
 In the eyepiece
- 
**SD card**  
 For data storage
- 
**USB 2.0 digital camera**  
 For direct transmitting of the picture to a PC
- 
**USB 3.0 digital camera**  
 For direct transmitting of the picture to a PC
- 
**WiFi data interface:**  
 For transmitting of the picture to a mobile display device
- 
**HDMI digital camera**  
 For direct transmitting of the picture to a display device
- 
**PC software**  
 To transfer the measurements from the device to a PC.
- 
**Automatic temperature compensation**  
 For measurements between 10 °C and 30 °C
- 
**Protection against dust and water splashes IPxx:**  
 The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013
- 
**Battery operation**  
 Ready for battery operation. The battery type is specified for each device.
- 
**Battery operation rechargeable**  
 Prepared for a rechargeable battery operation
- 
**Plug-in power supply**  
 230V/50Hz in standard version for EU. On request GB, AUS or USA version.
- 
**Integrated power supply unit**  
 Integrated in microscope. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.
- 
**Package shipment**  
 The time required to manufacture the product internally is shown in days in the pictogram.

## ABBREVIATIONS

- C-Mount** Adapter for the connection of a camera to a trinocular microscope
- FPS** Frames per second
- H(S)WF** High (Super) Wide Field (Eyepiece with high eye point for wearers of glasses)
- LWD** Long Working Distance
- N.A.** Numerical Aperture
- SLR camera** Single-Lens Reflex camera
- SWF** Super Wide Field (Field number at least  $\varnothing$  23 mm for 10 $\times$  eyepiece)
- W.D.** Working Distance
- WF** Wide Field (Field number up to  $\varnothing$  22 mm for 10 $\times$  eyepiece)