

Precision balances KERN PWS



**NEW** **STANDARD**  
★★



High-resolution precision balance made of stainless steel with IP protection

### Features

- High-quality tuning fork measuring system for steady weight values and continuous weighing
- Thanks to the stainless steel design of the housing and platform with smooth surface, the scale is rust-free and easy to clean
- Dust and spray protection to IP65 (in accordance with EN 60529)
- Level indicator and levelling feet for precise levelling of the scale, fitted as standard
- RS-232 data interface for connection to a printer as standard
- Protective working cover included with delivery


### Technical data

- Large LCD display, digit height 16,5 mm
- Dimensions weighing surface, stainless steel
  - A ø 140 mm
  - B WxD 190x190 mm, see larger picture
- Overall dimensions WxDxH 310x208x87 mm
- Net weight approx. 1,6 kg
- Permissible ambient temperature 10 °C/30 °C

### Accessories

- Loop for underfloor weighing, KERN EG-A07
- Software BalanceConnection, for flexible recording or transmission of measured values, in particular also to Microsoft® Excel or Access as well as transfer of this data to other Apps and programs, for more details see the internet, Scope of supplies: 1 CD, 1 license, KERN SCD-4.0
- Battery insert, enables mobile weighing, Operating time 200 h, KERN PWS-A01
- Interface cable RS-232 to connect an external device, Dust and spray protection to IP65 (in accordance with EN 60529), KERN PWS-A02
- Further details, plenty of further accessories and suitable printers see *Accessories*

#### STANDARD

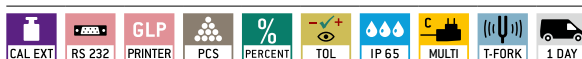

#### OPTION



Model	Weighing capacity [Max]	Readability [d]	Reproducibility	Linearity	Weighing plate	Option
KERN	g	g	g	g		DAkKS Calibr. Certificate
PWS 800-2	820	0,01	0,01	± 0,01	A	DAkKS KERN 963-127
PWS 3000-1	3200	0,1	0,1	± 0,1	B	963-127
PWS 8000-1	8200	0,1	0,1	± 0,1	B	963-128

EXPECTED TO BE AVAILABLE FROM MARCH 2023

- 
**Internal adjusting:**  
 Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)
- 
**Adjusting program CAL:**  
 For quick setting up of the balance's accuracy. External adjusting weight required
- 
**Easy Touch:**  
 Suitable for the connection, data transmission and control through PC or tablet.
- 
**Memory:**  
 Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.
- 
**Alibi memory:**  
 Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.
- 
**KERN Universal Port (KUP):**  
 allows the connection of external KUP interface adapters, e.g. RS-232, RS-485, SB, Bluetooth, WLAN, Analogue, Ethernet etc. for the exchange of data and control commands, without installation effort
- 
**Data interface RS-232:**  
 To connect the balance to a printer, PC or network
- 
**RS-485 data interface:**  
 To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible
- 
**USB data interface:**  
 To connect the balance to a printer, PC or other peripherals
- 
**Bluetooth\* data interface:**  
 To transfer data from the balance to a printer, PC or other peripherals
- 
**WiFi data interface:**  
 To transfer data from the balance to a printer, PC or other peripherals
- 
**Control outputs (optocoupler, digital I/O):**  
 To connect relays, signal lamps, valves, etc.
- 
**Analogue interface:**  
 to connect a suitable peripheral device for analogue processing of the measurements
- 
**Interface for second balance:**  
 For direct connection of a second balance
- 
**Network interface:**  
 For connecting the scale to an Ethernet network
- 
**KERN Communication Protocol (KCP):**  
 It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems
- 
**GLP/ISO log:**  
 The balance displays weight, date and time, independent of a printer connection
- 
**GLP/ISO log:**  
 With weight, date and time. Only with KERN printers.
- 
**Piece counting:**  
 Reference quantities selectable. Display can be switched from piece to weight
- 
**Recipe level A:**  
 The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out
- 
**Recipe level B:**  
 Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display
- 
**Totalising level A:**  
 The weights of similar items can be added together and the total can be printed out
- 
**Percentage determination:**  
 Determining the deviation in % from the target value (100 %)
- 
**Weighing units:**  
 Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details
- 
**Weighing with tolerance range:**  
 (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model
- 
**Hold function:**  
 (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value
- 
**Protection against dust and water splashes IPxx:**  
 The type of protection is shown in the pictogram.
- 
**Suspended weighing:**  
 Load support with hook on the underside of the balance
- 
**Battery operation:**  
 Ready for battery operation. The battery type is specified for each device
- 
**Rechargeable battery pack:**  
 Rechargeable set
- 
**Universal plug-in power supply:**  
 with universal input and optional input socket adapters for  
 A) EU, CH, GB  
 B) EU, CH, GB, USA  
 C) EU, CH, GB, USA, AUS
- 
**Plug-in power supply:**  
 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available
- 
**Integrated power supply unit:**  
 Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
- 
**Weighing principle: Strain gauges**  
 Electrical resistor on an elastic deforming body
- 
**Weighing principle: Tuning fork**  
 A resonating body is electromagnetically excited, causing it to oscillate
- 
**Weighing principle: Electromagnetic force compensation**  
 Coil inside a permanent magnet. For the most accurate weighings
- 
**Weighing principle: Single cell technology:**  
 Advanced version of the force compensation principle with the highest level of precision
- 
**Verification possible:**  
 The time required for verification is specified in the pictogram
- 
**DAkkS calibration possible (DKD):**  
 The time required for DAkkS calibration is shown in days in the pictogram
- 
**Factory calibration (ISO):**  
 The time required for Factory calibration is shown in days in the pictogram
- 
**Package shipment:**  
 The time required for internal shipping preparations is shown in days in the pictogram
- 
**Pallet shipment:**  
 The time required for internal shipping preparations is shown in days in the pictogram

\*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.