

IP protected bench scale KERN WTB



Ideal for food industry, thanks to its IP65 protection and special checkweighing display, now also with EC type approval [M]

Features

- Ideal for the increased hygienic requirements in the food industries
- Your support in a HACCP-compliant quality system
- **1** Protection against dust and water splashes IP65 (only when using rechargeable battery pack)
- Very fast display: steady weight values within 2 s
- **2** Second display on the back of the balance
- Compact size, practical for small spaces
- **3** Innovative weighing with tolerance range (Checkweighing): The colour of the display changes depending on the weight (too light/ok/too heavy) and is an aid to help with portioning, dosing and grading
- High mobility: thanks to rechargeable battery operation (optional), compact,

lightweight construction, it is suitable

- for the use in several locations
- Recessed grips on the underside of the scale for easy transportation
- Protective working cover included with delivery


Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions weighing surface, stainless steel, WxD 263x203 mm
- Overall dimensions WxDxH 286x316x126,5 mm
- Net weight approx. 3,2 kg
- Permissible ambient temperature -10 °C/40 °C


Accessories

- Protective working cover, scope of delivery 5 items, KERN RFB-A01S05
- Internal rechargeable battery pack, operating time up to 50 h without backlight, charging time approx. 12 h, KERN WTB-A01N
- **4** Tare pan made from stainless steel, overall dimensions WxDxH, 400x300x45 mm, KERN RFS-A02

STANDARD

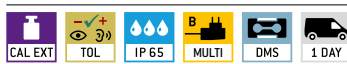

OPTION

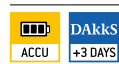

FACTORY



Model	Weighing capacity [Max] kg	Readability [d] g	Verification value [e] g	Minimal load [Min] g	Option	
					Verification M KERN	DAkkS Calibr. Certificate DAkkS KERN
KERN						
WTB 1K-4N	1,5	0,2	-	-	-	963-127
WTB 3K-4N	3	0,5	-	-	-	963-127
WTB 6K-3N	6	1	-	-	-	963-128
WTB 10K-3N	15	2	-	-	-	963-128
WTB 30K-3N	30	5	-	-	-	963-128

Note: For applications that require verification, please order verification at the same time, initial verification at a later date is not possible.
Verification at the factory, we need to know the full address of the location of use.

WTB 1K-4NM	1,5	0,5	0,5	10	965-227	963-127
WTB 3K-3NM	3	1	1	20	965-227	963-127
WTB 6K-3NM	6	2	2	40	965-228	963-128
WTB 10K-3NM	15	5	5	100	965-228	963-128
WTB 30K-2NM	30	10	10	200	965-228	963-128

- 
Internal adjusting:
 Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)
- 
Adjusting program CAL:
 For quick setting up of the balance's accuracy. External adjusting weight required
- 
Easy Touch:
 Suitable for the connection, data transmission and control through PC or tablet.
- 
Memory:
 Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.
- 
Alibi memory:
 Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.
- 
KERN Universal Port (KUP):
 allows the connection of external KUP interface adapters, e.g. RS-232, RS-485, SB, Bluetooth, WLAN, Analogue, Ethernet etc. for the exchange of data and control commands, without installation effort
- 
Data interface RS-232:
 To connect the balance to a printer, PC or network
- 
RS-485 data interface:
 To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible
- 
USB data interface:
 To connect the balance to a printer, PC or other peripherals
- 
Bluetooth* data interface:
 To transfer data from the balance to a printer, PC or other peripherals
- 
WiFi data interface:
 To transfer data from the balance to a printer, PC or other peripherals
- 
Control outputs (optocoupler, digital I/O):
 To connect relays, signal lamps, valves, etc.
- 
Analogue interface:
 to connect a suitable peripheral device for analogue processing of the measurements
- 
Interface for second balance:
 For direct connection of a second balance
- 
Network interface:
 For connecting the scale to an Ethernet network
- 
KERN Communication Protocol (KCP):
 It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems
- 
GLP/ISO log:
 The balance displays weight, date and time, independent of a printer connection
- 
GLP/ISO log:
 With weight, date and time. Only with KERN printers.
- 
Piece counting:
 Reference quantities selectable. Display can be switched from piece to weight
- 
Recipe level A:
 The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out
- 
Recipe level B:
 Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display
- 
Totalising level A:
 The weights of similar items can be added together and the total can be printed out
- 
Percentage determination:
 Determining the deviation in % from the target value (100 %)
- 
Weighing units:
 Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details
- 
Weighing with tolerance range:
 (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model
- 
Hold function:
 (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value
- 
Protection against dust and water splashes IPxx:
 The type of protection is shown in the pictogram.
- 
Suspended weighing:
 Load support with hook on the underside of the balance
- 
Battery operation:
 Ready for battery operation. The battery type is specified for each device
- 
Rechargeable battery pack:
 Rechargeable set
- 
Universal plug-in power supply:
 with universal input and optional input socket adapters for
 A) EU, CH, GB
 B) EU, CH, GB, USA
 C) EU, CH, GB, USA, AUS
- 
Plug-in power supply:
 230V/50Hz in standard version for EU, CH. On request GB, USA or AUS version available
- 
Integrated power supply unit:
 Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
- 
Weighing principle: Strain gauges
 Electrical resistor on an elastic deforming body
- 
Weighing principle: Tuning fork
 A resonating body is electromagnetically excited, causing it to oscillate
- 
Weighing principle: Electromagnetic force compensation
 Coil inside a permanent magnet. For the most accurate weighings
- 
Weighing principle: Single cell technology:
 Advanced version of the force compensation principle with the highest level of precision
- 
Verification possible:
 The time required for verification is specified in the pictogram
- 
DAkkS calibration possible (DKD):
 The time required for DAkkS calibration is shown in days in the pictogram
- 
Factory calibration (ISO):
 The time required for Factory calibration is shown in days in the pictogram
- 
Package shipment:
 The time required for internal shipping preparations is shown in days in the pictogram
- 
Pallet shipment:
 The time required for internal shipping preparations is shown in days in the pictogram

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.