



6 Block Dry Block Heaters



Process More Samples with Large Capacity Block Heaters that Deliver Precise, Reliable Results

Versatile Dry Block Heaters are ideal for applications that require temperature stability. All units hold one optional interchangeable modular block with over 40 choices. Blocks are easy to use, no tools required. Close tube-to-block contact enables maximum heat retention, resulting in efficient heating. Digital models offer exceptional temperature uniformity and stability for repeatable results, while non-digital models are an economical choice.

Unique Features:

- Unique design for high throughput applications. The six-block capacity allows you the most flexibility to process multiple samples simultaneously with one device.
- Ensure temperature accuracy. Digital models include a unique single point calibration feature to match the block heater's display to a known external temperature standard for 3 temperature setpoints.
- Digital models feature easy to use touchpad controls with independent displays for temperature and time. Robust non-digital models are controlled with dual temperature adjustment knobs.

6 Block Dry Block Heaters

Model	HB6AL	HB6DG
Control	Non-Digital	Digital
Temperature Range	Low Ambient + 5 °C – 100 °C; High 75 °C – 150 °C	Ambient + 5 °C – 120 °C
Stability @ 37 °C	+/-2 °C (120V) +/-2.5 °C (230V)	+/-0.1 °C (120V) +/-0.2 °C (230V)
Uniformity within the block @ 37 °C	+/-0.1 °C (120V) +/-0.4 °C (230V)	+/-0.1 °C (120V) +/-0.2 °C (230V)
Heat Up Time to 100 °C	≤ 30 min	
Temperature Calibration	NA	.
Number of Blocks	6	
Overshoot Protection	NA	10 °C above set point
Working Environment	18 °C – 33 °C, 20% – 80% relative humidity, non-condensing	
Power	230V, 1.79A, 50/60Hz	
Power Consumption	410 W	
Dimensions (W x L x H)	203 x 531 x 89 mm	
Net Weight	4.4 kg	

Other Standard Features and Equipment

Touchpad control and audible alarm (digital model), painted steel housing, digital model includes support rod holder with locking knob, detachable 3-wire cord and plug (included)

Compliance


- **Product Safety:** EN 61010-1, EN 61010-2-051
- **Electromagnetic Compatibility:** EN 61326-1 Class A, Industrial Environments
- **Compliance Marks:** CE; RCM; UKCA; TÜV SÜD

Accessories

Module Block 10 mm Test Tube, 24 Wells	30400151	Sand Bath, 6 Block Unit	30400176
Module Block 12/13 mm, 20 Wells	30400152	Sand, 0.5 kg Bag	30400177
Module Block 15/16 mm, 12 Wells	30400153	Stainless Steel Shot, 0.5 kg	30400178
Module Block 20 mm Test Tube, 8 Wells	30400154	Cover Low-Temperature 6 Block	30400181
Module Block 25 mm Test Tube, 6 Wells	30400155	Module Block For Vials 12 mm, 20 Wells	30400182
Module Block Combo - Test Tube Combination Block:		Module Block For Vials 15 mm, 20 Wells	30400183
6 mm, 12/13 mm, 25 mm	30400156	Module Block For Vials 17 mm, 12 Wells	30400184
Module Block 0.5 mL Microtaper, 30 Wells	30400157	Module Block For Vials 19 mm, 12 Wells	30400185
Module Block 6 mm Test Tube, 30 Wells	30400158	Module Block For Vials 21 mm, 9 Wells	30400186
Module Block 1.5 mL Microtaper, 20 Wells	30400159	Module Block For Vials 23 mm, 8 Wells	30400187
Module Block Solid	30400160	Module Block For Vials 25 mm, 8 Wells	30400188
Module Block Cuvette - 12 Cuvet, 2 Wells	30400161	Module Block For Vials 28 mm, 6 Wells	30400189
Module Block 11.5/1.5 mL Microcentrifuge, 20 Wells		Module Block For Vials 16 mm, 15 Wells	30400190
.	30400162	Module Block, 2.0 mL Tube, 20 Wells	30400191
Module Block 250 Microliter/6 mm, 30 Wells	30400163	Module Block, 2.0 mL Corning Tube, 20 Wells	30400192
Module Block Titer Plate	30400164	Module Block Combo - Centrifuge Tube Configuration:	
Module Block 12/13 mm, 16 Wells	30400165	50 mL-2 Wells, 15 mL-3 Wells, 1.5 mL-4 Wells	30400193
Module Block Solid Double	30400166	Module Block Combo - Micro Tube Combination:	
Module Block 35 mm Test Tube, 4 Wells	30400167	2 mL-5 Wells, 1.5 mL-10 Wells, 0.5 mL-6 Wells	30400194
Module Block 50 mL Conical Centrifuge, 5 Wells	30400168	Module Block 17/18 mm, 12 Wells	30400195
Module Block 0.2 mL Micro, 80 Wells	30400169	Stainless Steel Probe Kit*	30400246
Module Block 0.2 mL Micro, 64 Wells	30400170		
Module Block 0.2 mL 96 Hole Double Block	30400171		

*NA for HB6AL

Dimensions

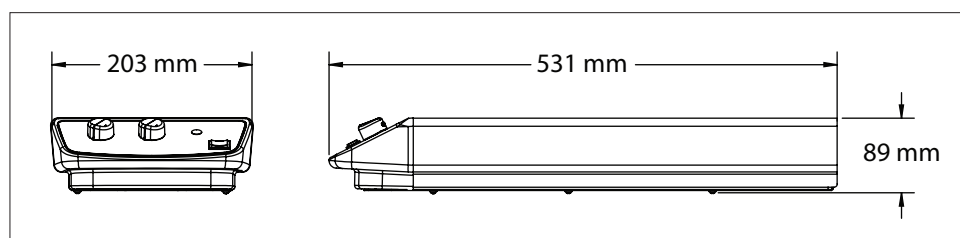

775641_B 20221130 © Copyright OHAUS Corporation

